

Instructions for use

miniRITE R

Oticon Opn S™

Oticon Opn Play™

Oticon Ruby



Made for
iPhone | iPad | iPod

oticon
life-changing technology

Model overview

This booklet is valid for the Oticon Opn S™, Oticon Opn Play™ and Oticon Ruby families in the following hearing aid models:

FW 8

- | | |
|---|---------------------------|
| <input type="checkbox"/> Oticon Opn S 1 | GTIN: (01) 05707131374967 |
| <input type="checkbox"/> Oticon Opn S 2 | GTIN: (01) 05707131374974 |
| <input type="checkbox"/> Oticon Opn S 3 | GTIN: (01) 05707131374981 |

- | | |
|--|---------------------------|
| <input type="checkbox"/> Oticon Opn Play 1 | GTIN: (01) 05707131374998 |
| <input type="checkbox"/> Oticon Opn Play 2 | GTIN: (01) 05707131375001 |

FW 1

- | | |
|--|---------------------------|
| <input type="checkbox"/> Oticon Ruby 1 | GTIN: (01) 05707131378491 |
| <input type="checkbox"/> Oticon Ruby 2 | GTIN: (01) 05707131378507 |

The following speakers are available for the above models:

- Speaker 60
- Speaker 85
- Speaker 100 (Power Instrument)

- Power flex mould speaker 100 (Power Instrument)
- Power flex mould speaker 105 (Power Instrument)

Introduction to this booklet

This booklet gives you guidance on how to use and maintain your new hearing aid. Please read the booklet carefully, including the **Warnings** section. This will help you to get the most out of your new hearing aid.

Your hearing care professional has adjusted the hearing aid to meet your needs. If you have further questions, please contact your hearing care professional.

For your convenience, this booklet contains a navigation bar to help you easily navigate through the different sections.

| **About** | [Start-up](#) | [Handling](#) | [Options](#) | [Tinnitus](#) | [Warnings](#) | [More info](#) |

Intended use

Intended use	The hearing aid is intended to amplify and transmit sound to the ear.
Indications for use	Impaired hearing within mild to severe-to-profound hearing loss.
Intended user	Adults and children older than 36 months.
User environment	Indoor and outdoor.
Contraindications	Users of active implants must pay special attention when using the hearing aid. For more information read the Warnings section.
Clinical benefits	The hearing aid is designed to provide better speech understanding to help ease communication with the aim of improving quality of life.

IMPORTANT NOTICE

The hearing aid amplification is uniquely adjusted and optimised to your personal hearing capabilities during the hearing aid fitting performed by your hearing care professional.

In short

Charging

Ensure you charge the hearing aid every night.

The LED on the hearing aid will turn RED/ORANGE when the hearing aid is charging and GREEN when fully charged.

It takes 3 hours to fully charge the hearing aid.

ON/OFF

The hearing aid turns ON when you remove it from the charger.

The hearing aid turns OFF when you place it in the charger.

Always leave the charger connected to a power source when the hearing aid is seated in the charging port. Switching off the charger will make the hearing aid turn on and start using battery.

Table of contents

About

Your hearing aid, speaker and earpiece	10
Hearing aid	12

Start-up

Charging time	14
Battery performance	15
Turn hearing aid ON and OFF using charger	16
Turn hearing aid ON and OFF using push button	17
Battery low indication	18
Identify left and right hearing aid	19
Put on hearing aid	20

Continues on next page

Table of contents

Handling

Cleaning	22
Replace standard earpieces	24
ProWax miniFit filter	26
Clean customised earpieces	28
Replace ProWax filter	29
Storage of the hearing aid	30
Flight mode	32

Options

Optional features and accessories	34
Change volume	35
Change program	36
Mute	37
Use hearing aid with iPhone and iPad	38
Wireless accessories and other options	40

Tinnitus

Tinnitus SoundSupport™ (optional)	42
Sound and volume adjustments	44
Limitation on use time	47
Tinnitus SoundSupport warnings	50


Warnings

General warnings	52
------------------	----

More info



Troubleshooting	60
Water & dust resistant (IP68)	64
Conditions of use	65
Technical information	66
Australia only	76
Warranty	77
Your individual hearing aid settings	78
Sound and LED indicators	80

Your hearing aid, speaker and earpiece

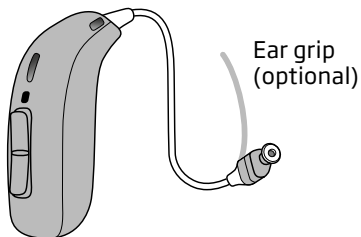


The hearing aid uses one of the following speakers:

**Standard
speakers**

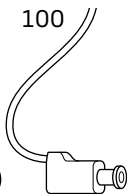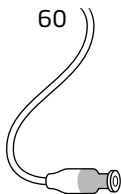

**Power flex
mould
speakers**

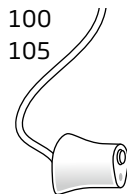

The speakers uses one of the following earpieces:

Standard earpieces



Open dome



Bass dome, single vent



Bass dome, double vent



Power dome



Grip Tip

Available in small and large, left and right, with or without vent.

Customised earpieces



Micro mould



LiteTip

Dome sizes




* only as open dome for speaker 60

** not for open dome

For details on replacing the dome, see the "Replace standard earpieces" section.

Hearing aid

What it is and does

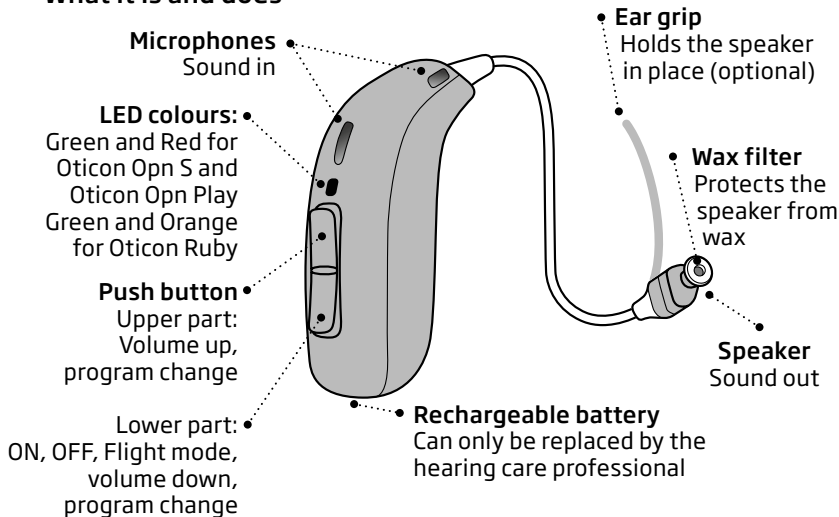

Charging time

Ensure you fully charge your hearing aid before first time use.

Charge your hearing aid every night. That ensures you start your day with a fully charged hearing aid.

If your hearing aid is completely drained the normal charging time is:

3 hours	1 hour	0.5 hour
Fully charged	50% charged	25 % charged

The charging process will automatically stop when the battery is fully charged. Leave the hearing aid in the charger when you do not use it. Remember to keep the charger connected to a power source.

Charging time may vary depending on the remaining capacity of the battery and between the left and right hearing aid.

For how to use your charger see the charger **Instructions** for use.

Battery performance

The daily battery performance varies depending on your individual use and hearing aid settings.

Streaming sound from e.g. TV, mobile phone or ConnectClip can influence daily battery performance.

If your hearing aid runs out of battery, place it in the charger for a re-charge. Additional usage time cannot be achieved by trying to restart the hearing aid.

For best charging condition, the room temperature should be between +10°C to +35°C.

If you experience that your hearing aid does not perform a full day, it might be time to have the battery replaced. Contact your hearing care professional.

Turn hearing aid **ON** and **OFF** using charger

Your hearing aid automatically turns **ON** when removed from the charger.

The hearing aid LED turns **GREEN** after 6 seconds, confirming that it is ready for use. You may hear a start-up jingle.


Your hearing aid automatically turns **OFF** when placed in the charger and starts charging. The LED turns **RED/ORANGE**.

IMPORTANT NOTICE

Always leave the charger connected to a power source when the hearing aid is seated in the charging port. Switching off the charger will make the hearing aid turn on and start using battery.

Turn hearing aid ON and OFF using push button

The hearing aid can be turned ON and OFF using the push button.

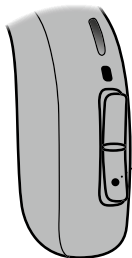

ON

Press and hold the lower part of the push button for approximately 2 seconds.

Release the button and wait until the hearing aid LED turns GREEN.

The hearing aid is now turned ON. You may hear a short start-up jingle.

OFF

Press and hold the lower part of the button for 3 seconds until the hearing aid LED turns RED/ORANGE.

The hearing aid plays 4 descending tones.

Release the push button and the hearing aid is turned OFF.

See overview "Sound and LED indications".

Battery low indication

When the battery is running low, you will hear three short beeps. This will leave you approximately two hours before the hearing aid stops working. The beeps will be repeated every 30 minutes. Just before the battery runs out you will hear four descending tones.

If streaming audio to your hearing aid (from e.g. TV or mobile phone) you have approximately one hour before the hearing aid stops working, when you hear the three short beeps. Stopping streaming will prolong the battery performance.



Three alternate tones
= The battery is running low



Four descending tones
= The battery has run out


Optional LED

Continuous red/orange flashes indicate low battery.

Identify left and right hearing aid

It is important to distinguish between the left and the right hearing aid, as they may be programmed differently.

You can find left/right colour indicators on the hearing aid itself and on 60 and 85 speakers as shown. Indicators (either L or R) can also be found on 100 speakers and some earpieces.

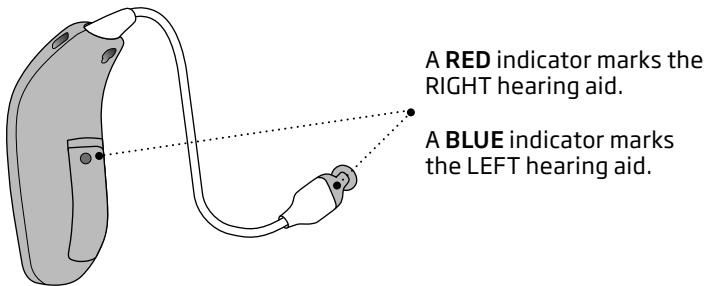

Put on hearing aid


Step 1



Place the hearing aid behind your ear.


The speaker should always be used with an earpiece attached. Use only parts designed for your hearing aid.

Step 2

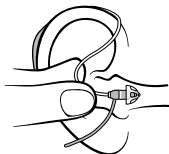

Hold the bend of the speaker wire between your thumb and index finger. The earpiece should point toward the opening of the ear canal.

Step 3

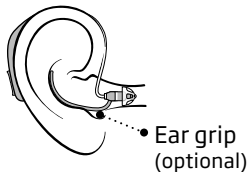


Gently push the earpiece into your ear canal until the speaker wire sits close to your head.

If the speaker has an ear grip, place it in the ear so it follows the contour of the ear.

Cleaning

The MultiTool contains a brush and wire loop for cleaning and removing earwax. If you need a new MultiTool, please contact your hearing care professional.

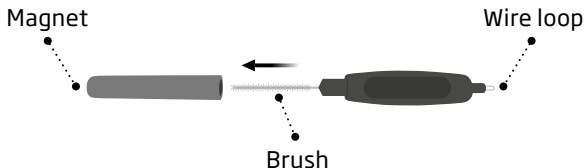

IMPORTANT NOTICE


The MultiTool has a built-in magnet. Keep the MultiTool 30 cm away from credit cards and other magnetically sensitive devices.

When handling the hearing aid, hold it over a soft surface to avoid damage if you drop it.

Clean the microphone openings

Use the brush of the MultiTool to carefully brush debris away from the openings. Carefully brush the surface around the opening.

Make sure that no parts of the MultiTool are squeezed into the microphone openings by force. This may damage the hearing aid.

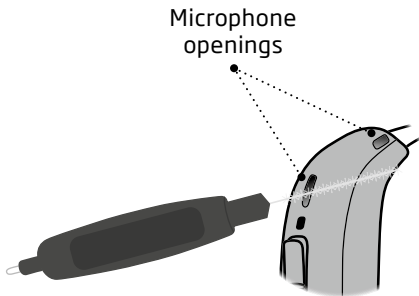


IMPORTANT NOTICE

Use a soft, dry cloth to clean the hearing aid. The hearing aid must never be washed or immersed in water or other liquids.

Replace standard earpieces


The standard earpiece (dome and Grip Tip) should not be cleaned. If the earpiece is filled with wax, replace it with a new one. Grip Tip needs to be replaced at least once a month.

Step 1

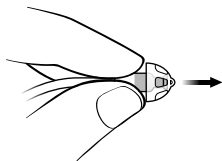


Hold on to the speaker and pull off the earpiece.

Step 2

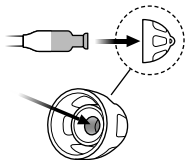

Insert the speaker exactly into the middle of the earpiece to obtain a secure attachment.

Step 3

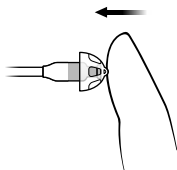

Push firmly to ensure that the earpiece is fastened securely.


IMPORTANT NOTICE

If the earpiece is not on the speaker when removed from the ear, the earpiece may still be in the ear canal. Consult your hearing care professional for further instructions.

ProWax miniFit filter

The speaker has a white wax filter attached to the end where the earpiece is attached. The wax filter keeps wax and debris from damaging the speaker. Replace the filter when clogged, or if the hearing aid does not sound normal. Alternatively, contact your hearing care professional.

Remove the earpiece from the speaker before replacing the wax filter.

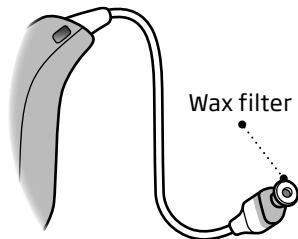


IMPORTANT NOTICE

Always use the same type of wax filter as was originally supplied with the hearing aid. If you are in any doubt about the use or replacement of wax filters, contact your hearing care professional.


Replace ProWax miniFit filter

1. Tool

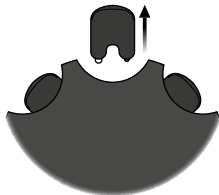


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new wax filter.

2. Remove

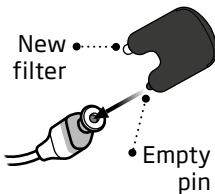

Push the empty pin into the wax filter in the speaker and pull it out.

3. Insert

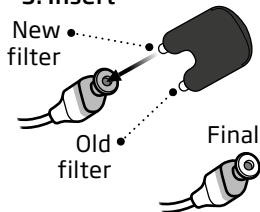

Insert the new wax filter using the other pin, remove the tool, and throw it out.

Note:


If you use a mould or LiteTip, your hearing care professional must replace the wax filter in the speaker.

Clean customised earpieces

The earpiece should be cleaned regularly.


The earpiece has a white wax filter. The filter keeps wax and debris from damaging the speaker.

Replace the filter when clogged, or if the hearing aid does not sound normal. Alternatively, contact your hearing care professional.

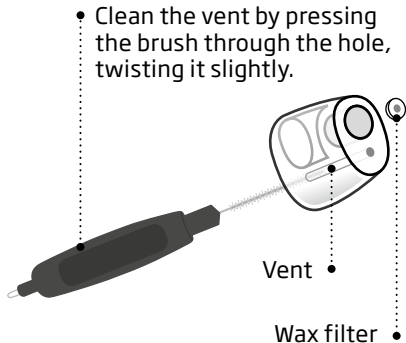


Replace ProWax filter

1. Tool

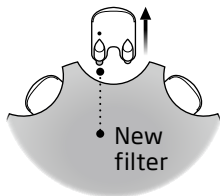


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new wax filter.

2. Remove

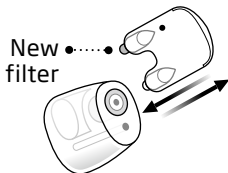

Push the empty pin into the wax filter in the earpiece and pull it out.

3. Insert

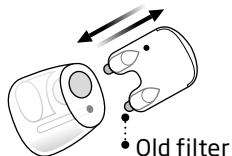

Insert the new wax filter using the other pin, remove the tool, and throw it out.

Storage of the hearing aid

The charger is the best storage for the hearing aid. Anytime your hearing aid is not in use, place it in the charger for charging. This will ensure your hearing aid is always charged.

To ensure the longest life of the rechargeable battery in the hearing aid do not expose to excessive heat. For example, do not leave the hearing aid in the sun in a window or in a car, even if the hearing aid is in the charger.

IMPORTANT NOTICE

Always leave the charger connected to a power source when the hearing aid is seated in the charging port. Switching off the charger will make the hearing aid turn on and start using power.

Long term storage (weeks, months)


Before storing the hearing aid for a longer period (more than 14 days) ensure you fully charge the hearing aid. Then turn the hearing aid off. This way the battery can be charged again.

It is necessary to fully charge the hearing aid every 6 months to protect the rechargeable battery. If a stored hearing aid is not charged within a 6 months time frame, the rechargeable battery will need to be replaced.

To see how to place the hearing aid in the charger see the Charger Instructions for use.

Flight mode

To activate and deactivate flight mode press and hold the push button (lower part) for 7 seconds.

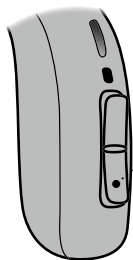

Beep: The hearing aid plays 4 descending beeps followed by 4 beeps. This indicates that flight mode is activated or deactivated.

LED flight mode activated: A long red/orange blink followed by green, red/orange, red/orange. The LED will blink green, red/orange, red/orange a few times. This confirms that flight mode is activated.

LED flight mode deactivated: A long red/orange blink followed by two long green blinks. This confirms that flight mode is deactivated.

See overview "Sound and LED indications".

Pressing the button on one hearing aid will activate flight mode on both hearing aids.

When flight mode is activated, Bluetooth® is turned off. The hearing aid will still be working.

Optional features and accessories


The features and accessories described in the following pages are optional. Please contact your hearing care professional for more information about features and accessories.

If you experience difficult listening situations, a special program may be helpful. These are programmed by your hearing care professional.

Write down hearing situations in which you may need help.

Change volume


The push button allows you to adjust the volume.
You will hear a beep when you increase or decrease the volume.

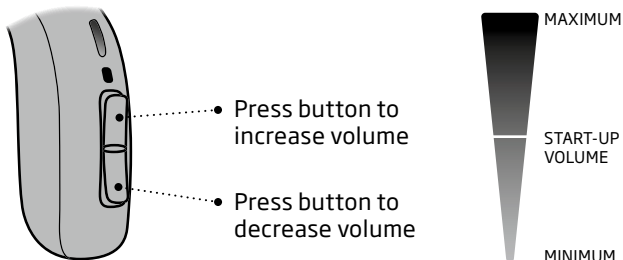

Change program

Your hearing aid can have up to four different programs. These are programmed by your hearing care professional. You will hear one to four tones when you change program depending on the program.

See the Sound and LED indicator section.

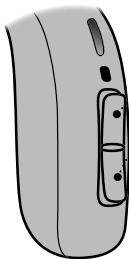

Press the push button to switch between programs.

The program cycle switches one program forward, for example program 1 to 2 or program 4 to 1.

Mute

Use the mute function if you need to silence the hearing aid. Your hearing aid can be muted by using one of the following optional devices:

Oticon ON App

ConnectClip

Remote Control 3.0

How to unmute your hearing aid

The hearing aid can be unmuted using one of the optional devices or by applying a short press to the upper or lower part of the button on the hearing aid.

IMPORTANT NOTICE

Do not use the mute function as an off switch, as the hearing aid still draws current from the battery in this mode.

Use hearing aid with iPhone and iPad

Your hearing aid is Made for iPhone® and allows for direct communication and control with compatible iPhone, iPad® or iPod touch® devices. For assistance in using these products with your hearing aid, please contact your hearing care professional.



To see how to pair your hearing aid with your iPhone please visit www.oticon.com.au/pairing

For information on compatibility, please visit www.oticon.com.au/compatibility

Use of the Made for Apple badge means that an accessory has been designed to connect specifically to the Apple products identified in the badge, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that use of this Accessory with iPod, iPhone, or iPad may affect wireless performance.

Wireless accessories and other options

As an enhancement to your wireless hearing aid, a range of accessories is available. These can enable you to hear and communicate better in many everyday situations.

ConnectClip

When ConnectClip is paired with your mobile phone, you can use the hearing aid as a hands-free headset. ConnectClip can also be used as a remote microphone

TV Adapter 3.0

TV Adapter streams sound directly from a TV or electronic audio device to your hearing aid.

Remote Control 3.0

Offers the ability to change program, adjust volume, or mute your hearing aid.

Oticon ON App

Intuitive and discreet control of your hearing aid. For iPhone, iPad, iPod touch, and Android™ devices.

Telecoil program

Telecoil helps you hear better when using a telephone with a built-in loop or when you are in buildings with teleloop systems such as theatres, churches, or lecture rooms. This symbol or a similar sign is shown wherever a teleloop has been installed.



CROS - optional for Oticon Opn S 1, Opn S 2, Oticon Opn Play 1, Oticon Ruby 1 and Ruby 2

Solution for people with an unaidable hearing loss in one ear. Oticon CROS on the poorer ear transmits sound to the hearing aid on the better ear.

For more information visit www.oticon.com.au or contact your hearing care professional.

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, Google Play, and the Google Play logo are trademarks of Google LLC.

Tinnitus SoundSupport™ (optional)

Intended use of Tinnitus SoundSupport

Tinnitus SoundSupport is a tool intended to generate sounds to provide temporary relief for patients suffering from tinnitus as part of a tinnitus management program.

The target population is the adults (over 18 years old).

Tinnitus SoundSupport is targeted to licensed hearing care professionals (audiologists, hearing aid specialists, or otolaryngologists) who are familiar with the evaluation and treatment of tinnitus and hearing loss. Fitting of Tinnitus SoundSupport must be done by a hearing care professional participating in a tinnitus management program.

Guidelines for tinnitus sound generator users

These instructions contain information about Tinnitus SoundSupport, which may have been enabled in your hearing aids by your hearing care professional.

Tinnitus SoundSupport is a tinnitus management device intended to generate sound of sufficient intensity and bandwidth to help manage tinnitus.

Your hearing care professional will also be able to offer the appropriate follow-up care. It is important to follow his/her advice and directions regarding such care.

Prescription use only

Good health practice requires that the person reporting tinnitus has a medical evaluation by a licensed ear physician before using a sound generator. The purpose of such an evaluation is to ensure that any medically treatable condition that may cause tinnitus is identified and treated prior to using a sound generator.

Sound and volume adjustments

Tinnitus SoundSupport is programmed by your hearing care professional to match your hearing loss and preferences for tinnitus relief. It offers a number of different sound options. Together with your hearing care professional, you can select your preferred sound(s).

Tinnitus SoundSupport programs

Together with your hearing care professional you decide for which programs you may want to have Tinnitus SoundSupport activated. The sound generator can be activated in up to four different programs.

Volume adjustments with Tinnitus SoundSupport

When you select a hearing aid program for which Tinnitus SoundSupport is activated, your hearing care professional can only set the push button on your hearing aid to work as a volume control for the tinnitus relief sound.

Your hearing care professional will set the volume control for the sound generator in one of two ways:

- A) Change volume in each ear separately, or
- B) Change volume in both ears simultaneously.

miniRITE R

A) How to change Tinnitus SoundSupport volume in each ear separately

To **increase** volume (on one hearing aid only), use a short press on the upper part of the push button repeatedly until desired level is reached. To **decrease** volume (on one hearing aid only), use a short press on the lower part of the push button repeatedly until desired level is reached.

B) How to change Tinnitus SoundSupport volume in both ears simultaneously

You can use one hearing aid to increase/decrease the sound in both hearing aids. When changing the volume in one hearing aid, the volume on the other hearing aid will follow.

To **increase** volume, use a short press on the upper part of the push button repeatedly.

To **decrease** volume, use a short press on the lower part of the push button repeatedly.

To be filled out by your hearing care professional.

Limitation on use time

Daily use

The volume levels of Tinnitus SoundSupport can be set to a level which could lead to permanent hearing damage when used for a prolonged period of time. Your hearing care professional will advise you of the maximum amount of time per day you should use Tinnitus SoundSupport. It should never be used at uncomfortable levels.

See table “Tinnitus SoundSupport: Limitation on use” in section “Your individual hearing aid settings” at the end of this booklet to learn how many hours per day you can safely use the relief sound in your hearing aid.

Important information for hearing care professionals about Tinnitus SoundSupport

Device description

Tinnitus SoundSupport is a module function that can be enabled in the hearing aids by the hearing care professional.

Maximum wearing time

The wearing time of Tinnitus SoundSupport will decrease as you increase the level above 80 dB(A) SPL. The fitting software will automatically display a warning when the hearing aid exceeds 80 dB(A) SPL. See “Max wearing time indicator” next to the tinnitus fitting graph in the fitting software.

The volume control is deactivated

By default the volume control for the sound generator is deactivated in the hearing aid. Risk of noise exposure increases when the volume control is activated.

If the volume control is activated

A warning may be displayed if you activate the tinnitus volume control in the “Buttons & Indicators” screen. This occurs if the relief sound can be listened to at levels that may cause hearing damage. The “Max wearing time” table in the fitting software displays the number of hours the patient can safely use Tinnitus SoundSupport.

- The max wearing time for each program for which Tinnitus SoundSupport is activated.
- Write those values in the table: “Tinnitus SoundSupport: Limitation on use”, in the back of this booklet.
- Instruct your patient accordingly.

Tinnitus SoundSupport warnings

If your hearing care professional has activated the sound generator Tinnitus SoundSupport, please pay attention to the following warnings.

There are some potential concerns associated with the use of any sound generated by a tinnitus management device. Among them are the potential worsening of tinnitus, and/or a possible change in hearing thresholds.

Should you experience or notice a change in hearing or tinnitus, or any dizziness, nausea, headaches, heart palpitations, or possible skin irritation at the point of contact with the device, you should immediately discontinue use of the device and consult a medical, audiology, or other hearing care professional.

As with any device, misuse of the sound generator feature may cause potentially harmful effects. Care should be taken to prevent unauthorised use and to keep the device out of reach of children and pets.

Maximum wearing time

Always follow the maximum wearing time per day of the Tinnitus SoundSupport advised by your hearing care professional. Prolonged use may lead to worsening of your tinnitus or of your hearing loss.

General warnings

For your personal safety and to ensure correct usage, you should familiarise yourself fully with the following general warnings before using your hearing aid.

Consult your hearing care professional if you experience unexpected operations or serious incidents with your hearing aid during use or because of its use. Your hearing care professional will support you with issue handling and, if relevant, reporting to the manufacturer and/or the national authorities.

Please note that a hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions. Hearing aid is only a part of hearing habilitation and may need to be supplemented by auditory training and instruction in lipreading. Furthermore, note that in most cases, infrequent use

of a hearing aid does not permit a user to attain full benefit from it.

This hearing aid is supported by a non-removable rechargeable lithium-ion battery cell. Please ensure to charge the hearing aid and familiarise yourself with the safety and handling information related to rechargeable hearing aids.

Do not try to get access to the battery inserted in the hearing instrument. The battery must only be replaced by your hearing care professional.

Only charge the hearing aid with an Oticon charger. Other chargers risk destroying the hearing aid and battery.

If a battery or hearing aid is swallowed, see a doctor immediately.

Usage of hearing aids

Hearing aids should be used only as directed and adjusted by your hearing care professional. Misuse can result in sudden and permanent hearing loss.

Never allow others to wear your hearing aid, as incorrect usage could cause permanent damage to their hearing.

Choking hazards & risk of swallowing small parts

Hearing aids and their parts should be kept out of reach of children and anyone who might swallow these items or otherwise cause injury to themselves.

Explosives

The hearing aid is safe to use under normal usage conditions. The hearing aid has not been tested for compliance with international standards concerning explosive environments.

Therefore, do not use the hearing aid in environments with danger of explosions e.g. mines, oxygen rich environments or areas where flammable anaesthetics are handled.

Fatality hazards and risk of swallowing Lithium-ion batteries or placing them in the ear or nose

Never swallow Lithium-ion batteries neither place them in the ear or the nose as this may lead to serious injury or death in as little as 2 hours. This can be due to chemical burn, which can cause permanent damage to the nose or ear or potentially lead to perforation of the inner organs. If a Lithium-ion battery is swallowed or placed in the ear or nose, seek emergency medical treatment immediately.

Continues on next page

General warnings

Rechargeable battery

Do not attempt to open the hearing aid, as it may damage the battery.

Never attempt to replace the battery. If battery replacement is needed, please return your device to your hearing care professional. The service guarantee is void if there are signs of tampering.

The safety of recharging batteries using a USB connector is determined by the external signal source. When connected to external equipment plugged into a wall outlet, this equipment must comply with IEC 62368-1 or equivalent safety standards.

Dysfunction

Be aware of the possibility that your hearing aid may stop working without notice. Keep this in mind when you depend on warning sounds (e.g. when you are in traffic). The hearing aids may stop functioning, for instance if the batteries have expired or if the tubing is blocked by moisture or earwax.

Active implants

The hearing aid has been thoroughly tested and characterised for human health according to international standards for human exposure (Specific Absorption Ratio - SAR), induced electromagnetic power and voltages into the human body.

The exposure values are well below internationally accepted safety limits for SAR, Induced electromagnetic power and voltages into the human body defined in the standards for human health and coexistence with active medical implants such as pacemakers and heart defibrillators.

If you have an active brain Implant, please contact the manufacturer of your implantable device for information about the risk of disturbance.

The Autophone or MultiTool (which has a built-in magnet) should be kept at least 30 cm away from the implant, e.g. do not carry it in your breast pocket.

In general, please follow the guidelines recommended by the manufacturers of implantable defibrillators and pacemakers on use with magnets.

Cochlear implants

If you are using a cochlear implant (CI) on one ear and a hearing aid on the other ear, make sure to always keep your CI at least at a 1cm distance from your hearing aid. The magnetic field from the CI sound processors, coils, and magnets may permanently damage the speaker unit in your hearing aid. Never place the devices close together on a table e.g. when cleaning or changing batteries. Do not carry the CI and the hearing aid together in the same box.

X-ray/CT/MR/PET scanning, electrotherapy and surgery

Remove your hearing aid before X-ray, CT/MR/PET scanning, electrotherapy, surgery, etc. as your hearing aid may be damaged when exposed to strong fields.

Continues on next page

General warnings

Heat and chemicals

The hearing aid must never be exposed to extreme heat, e.g. left inside a parked car in the sun.

The hearing aid must not be dried in microwave ovens or other ovens.

The chemicals in cosmetics, hairspray, perfume, aftershave lotion, sunscreen lotion, and insect repellent can damage the hearing aid. Always remove your hearing aid before applying such products and allow time to dry before use.

Power instrument

Special care should be exercised in selecting, fitting and using a hearing aid whose maximum sound pressure capability exceeds 132 dB SPL (IEC 60138-4/IEC 711), as there may be risk of impairing the remaining hearing of the hearing aid user.

For information on whether your hearing aid is a power instrument, see the model overview.

Possible side effects

Hearing aids and earpieces may cause an accelerated accumulation of earwax.

The non-allergenic materials used in hearing aids may in rare cases cause a skin irritation or other side effects.

Please seek consultation with a physician if these conditions occur.

Use on aircraft

Your hearing aid has Bluetooth. On board an aircraft, the hearing aid must be put into flight mode to deactivate Bluetooth, unless Bluetooth is permitted by the flight personnel.

Connection to external equipment

The safety of the use of hearing aids when connected to external equipment (with an auxiliary input cable and/or USB cable and/or directly), is determined by the external signal source. When connected to external equipment plugged into a wall outlet, this equipment must comply with IEC 62368-1 or equivalent safety standards.

Incompatible non-supplied accessories

Only use accessories, transducers or cables supplied by the manufacturer of this device. Incompatible accessories may result in reduced electromagnetic compatibility (EMC) of your device.

Modification of hearing aids is not allowed

Changes or modifications not expressly approved by the manufacturer will void the warranty of the equipment.

Interference

The hearing aid has been thoroughly tested for interference according to the most stringent international standards.

Electromagnetic interference may occur in the vicinity of equipment with the symbol to the right. Portable and mobile RF (radio frequency) communications equipment can affect the performance of your hearing aid. If your hearing aid is affected by electromagnetic interference, move away from the source to reduce the interference



WARNING

New and used button/coin batteries are hazardous and should be kept away from children. Button/coin batteries can cause severe or fatal injuries in 2 hours or less if swallowed or placed inside any part of the body. Seek medical attention immediately if this is suspected.

If it is suspected a button/coin battery has been swallowed or otherwise placed inside any part of the body, you should contact the

Australian Poisons Information Centre
on 13 11 26

New Zealand National Poisons Centre
on 0800 764 766
immediately for 24/7 fast, expert advice.

Troubleshooting

Troubleshooting for Charger, see Charger Instructions for use

Symptom	Possible causes
No sound	Hearing aid is out of power
	Dead battery
	Clogged earpieces (dome, Grip Tip, or mould)
Intermittent or reduced sound	Clogged sound outlet
	Moisture
	Hearing aid is out of power
Squealing noise	Hearing aid earpiece inserted incorrectly
	Earwax accumulated in ear canal
Beeping	If your hearing aid plays 8 beeps, 4 times consecutively, your hearing aid needs a microphone service check

If none of the above solutions work, consult your hearing care professional.

Solutions

Charge the hearing aid

Contact your hearing care professional

Clean mould

Replace wax filter, dome, or Grip Tip

Clean mould or replace wax filter, dome, or Grip Tip

Gently wipe the hearing aid and let it dry

Charge the hearing aid

Re-insert the earpiece

Have ear canal examined by your doctor

Contact your hearing care professional

Continues on next page

Troubleshooting

Troubleshooting for Charger, see Charger Instructions for use

Symptom	Possible causes
Hearing aid LED remains off when the hearing aid is placed in the charger	The charger is not turned on
	The hearing aid or charger is either too warm or too cold
	Charging is incomplete. The room temperature exceeds +35 °C, which prolongs the charging time. The charger has stopped charging to protect the battery
	The hearing aid is not correctly inserted into the charger
The hearing aid LED blinks red/orange when the hearing aid is placed in the charger	System error
Pairing issue with Apple device	Bluetooth connection failed
	Only one hearing aid paired

Solutions

Verify that the charger power supply is connected correctly

Move the charger and hearing aid to a location with a temperature between +5 °C to +40 °C

Reinsert the hearing aid into the charger. This will complete the charging within approximately 15 minutes

Check the charger slots for foreign objects

Contact your hearing care professional

- 1) Unpair your hearing aid.
- 2) Turn Bluetooth off and on again.
- 3) Turn off and turn back on the hearing aid.
- 4) Re-pair hearing aid (please visit www.oticon.com.au/pairing).

If none of the above solutions work, consult your hearing care professional.

Water & dust resistant (IP68)

Your hearing aid is dust tight and protected against ingress of water, which means it is designed to be worn in all daily life situations.

Therefore, you do not have to worry about sweat or getting wet in the rain.

Before charging the hearing aid make sure to wipe off any moisture.

Should your hearing aid come into contact with water and stop working, gently wipe off any water and let the hearing aid dry.

IMPORTANT NOTICE

Do not wear your hearing aid while showering or participating in water activities. Do not immerse your hearing aid in water or other liquids.

Conditions of use

Operational environment conditions	Temperature: +5°C to +40°C Humidity: 5% to 93%, non-condensing
Charging conditions	Temperature: +5°C to +40°C Humidity: 5% to 93%, non-condensing
Transportation and storage conditions	Temperature and humidity should not exceed the following limits for extended periods during transportation and storage: Transport: Temperature: -20°C to +60°C Humidity: 5% to 93%, non-condensing Storage: Temperature: -20°C to +30°C Humidity: 5% to 93%, non-condensing

Information about Charger: Condition of use, see Charger Instruction for use

Technical information

The hearing aid contains two radio technologies which are described below:

The hearing aid contains a radio transceiver using short range magnetic induction technology working at 3.84 MHz. The magnetic field strength of the transmitter is very weak and always below 15 nW (typically below $-40 \text{ dB}\mu\text{A/m}$ ($-12.20 \text{ dB}\mu\text{A/ft}$) at 10 meter (33 feet) distance).

The hearing aid also contains a radio transceiver using Bluetooth Low Energy (BLE) and a proprietary short-range radio technology, both working at 2.4 GHz. The 2.4 GHz radio

transmitter is weak and is always below 4 dBm e.i.r.p. in total radiated power.

The hearing aid complies with international standards concerning electromagnetic compatibility and human exposure. Only use your hearing aid in areas where wireless transmission is permitted.

Due to the limited space available on the hearing aid, relevant approval markings can be found in this booklet.

Additional information can be found in the "Technical data sheets" on www.oticon.com.au

The manufacturer declares that this hearing aid is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

Declaration of Conformity is available from the manufacturer.



Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark
www.oticon.global

CE 0543



Waste from electronic equipment must be handled according to local regulations.



 **Bluetooth®**



IP68

Australian Sponsor

Audmet Australia Pty Ltd
Suite 4, Level 4, Building B
11 Talavera Road
North Ryde NSW 2113

Description of symbols and abbreviations used in this booklet



Warnings

Text marked with a warning symbol must be read before using the device.



Manufacturer

The device is produced by the manufacturer whose name and address are stated next to the symbol. Indicates the medical device manufacturer, as defined in EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC.



CE mark

The device complies with Medical Device Directive 93/42/EEC. The four digit number indicates the identification of the notified body.



Electronic waste (WEEE)

Recycle hearing aids, accessories or batteries according to local regulations. Hearing aid users can also return electronic waste to their hearing care professional for disposal. Electronic equipment covered by Directive 2012/19/EU on waste and electrical equipment (WEEE).



Regulatory Compliance Mark (RCM)

The device complies with electrical safety, EMC and radio spectrum requirements for devices supplied to the Australian or New Zealand market.



Double insulation

Indicates that product power supply is double insulated. This means a very low risk of single fault failures i.e. single fault failures cannot occur.

**Indoor use**

Indicates that the product is qualified for indoor use.

IP68**IP code**

Indicates the class of protections against harmful ingress of water and particulate matter according to EN 60529:199. IP6X indicates total dust protection. IPX8 indicates the protection against the effects of continuous immersion in water.

**Bluetooth logo**

Registered trademark of Bluetooth SIG, Inc. where any use of such requires a license.

**Made for Apple badges**

Indicates that the device is compatible with iPhone, iPad and iPod touch.

**Hearing loop**

This logo incorporates the universal symbol for hearing assistance. The "T" signifies that a hearing loop is installed.

**Radio Frequency (RF) transmitter**

Your hearing aid contains an RF transmitter.

GTIN**Global Trade Item Number**

A globally unique 14-digit number used to identify medical device products including medical device software. GTIN in this booklet is related to medical device firmware. GTIN on regulatory packaging label is related to medical device hardware

Description of additional symbols used on labels



Keep dry

Indicates a medical device that needs to be protected from moisture.



Caution symbol

Consult instructions for use for warnings and cautions.



Catalogue number

Indicates the manufacturer's catalogue number so that the medical device can be identified.



Serial number

Indicates the manufacturer's serial number so that a specific medical device can be identified.



Medical Device

The device is a medical device.



Use-by date

Indicates the date after which the medical device is not to be used.



Li-ion

Battery recycling symbol

Li-ion Battery recycling symbol

**Temperature limit**

Indicates the temperature limits to which the medical device can be safely exposed

**Humidity limitation**

Indicates the range of humidity to which the medical device can be safely exposed.

**Button Battery Warning**

The device is intended to be used with button batteries. Refer to enclosed warnings before using the device.

International warranty

Your device is covered by an international limited warranty issued by the manufacturer for a period of 12 months from the date of delivery to you. This limited warranty covers manufacturing and material defects in the device itself, but not in accessories such as batteries, tubing, speakers, earpieces and filters, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, physical changes in your ear, damage due to foreign objects entering the device,

or incorrect adjustments are NOT covered by the limited warranty and may void it. Under this warranty, the manufacturer will choose, at its sole discretion, whether to repair the device, or replace it with an equivalent model. The above warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods. Your hearing care professional may have issued a warranty that goes beyond the clauses of this limited warranty. Please consult him/her for further information.

If you need service

Take your device to your hearing care professional, who may be able to sort out minor problems and make adjustments immediately. Your hearing care professional can assist you to obtain warranty service from the manufacturer. Your hearing care professional may charge a fee for their services.

Australia only

The rights given under this warranty are in addition to other rights and remedies under a law in relation to the good or services to which the warranty relates to. Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Australian customers should contact Oticon Australia at Suite 4, Level 4, Building B, 11 Talavera Road, North Ryde NSW 2113, telephone 1800 633 738 if you have any questions regarding this warranty.

Warranty

Certificate

Name of owner: _____

Hearing care professional: _____

Hearing care professional's address: _____

Hearing care professional's phone: _____

Purchase date: _____

Model left: _____ Serial no.: _____

Model right: _____ Serial no.: _____

Your individual hearing aid settings

To be filled out by your hearing care professional.

Tinnitus SoundSupport: Limitation on use			
<input type="checkbox"/>	No limitation on use		
	Program	Start-up volume (Tinnitus)	Max volume (Tinnitus)
<input type="checkbox"/>	1	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	2	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	3	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	4	Max _____ hours per day	Max _____ hours per day

Settings overview for your hearing aid

Left			Right	
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Change volume	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Change program	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Tinnitus SoundSupport	<input type="checkbox"/> Yes	<input type="checkbox"/> No
Volume control indicators				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at min/max volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps when changing volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at start-up volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
Battery indicators				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Low battery warning	<input type="checkbox"/> On	<input type="checkbox"/> Off

Sound and LED indicators









Different sounds and LED lights indicate the hearing aid status. The different indicators are listed on the following pages. For charger LED indicators, see Charger Instruction for use.

Your hearing care professional can set sound and LED indicators to match your preferences.

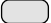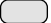Program	<input type="checkbox"/> Sound	<input type="checkbox"/> LED ¹⁾	When to use
1	1 tone	○	
2	2 tones	○ ○	
3	3 tones	○ ○ ○	
4	4 tones	○ ○ ○ ○	












○ Green flash

1) LED continuous or repeated three times with small pauses

ON/OFF	Sound	LED	LED comments
On	<input type="checkbox"/> Jingle	<input type="checkbox"/>  	Shown one time
Off	<input type="checkbox"/> 4 descending tones	<input type="checkbox"/> 	
Volume	Sound	LED	
Start-up volume	<input type="checkbox"/> 2 beeps	<input type="checkbox"/> 	
Minimum/maximum volume	<input type="checkbox"/> 3 beeps	<input type="checkbox"/> 	
Volume up/down	<input type="checkbox"/> 1 beep	<input type="checkbox"/> 	
Mute activated on ON App, ConnectClip or Remote Control 3.0		<input type="checkbox"/>  	Continuous or repeated three times


















 Long green flash
  Green flash
  Red/orange flash
  Long red/orange flash

Continues on next page

Accessories	<input type="checkbox"/> Sound	<input type="checkbox"/> LED	LED comments
TV Adapter	2 different tones	 	Continuous or repeated three times
ConnectClip remote microphone	2 different tones	 	
Flight mode	Sound	<input type="checkbox"/> LED	
Flight mode active	4 descending tones + short jingle	   	
Flight mode inactive	4 descending tones + short jingle	   ¹⁾	

 Long green flash  Green flash  Red/orange flash

1) Only available when three-time repetition is selected

Warnings	Sound	LED	LED comments
Low battery	 3 alternate tones	       	Continuously flashing
Battery shut down	4 descending tones		
Microphone service check needed	8 beeps repeated 4 times	    	Repeated four times
No light in the hearing aid LED when placed in the charger		Turned off	See Troubleshooting, solutions
The hearing aid LED flashes red/orange when the hearing aid is placed in the charger		  	Continuously flashing. See Troubleshooting, solutions

 Red/orange flash  Long red/orange flash



oticon
life-changing **technology**