

# Instructions for use

miniRITE R  
miniRITE T

Oticon Real™



Made for  
iPhone | iPad | iPod

Works with  
android

**oticon**  
life-changing technology

## Model overview

This booklet is valid for the following hearing aid models:

- miniRITE R
- miniRITE T

### FW 1.0

- |   |                           |
|---|---------------------------|
| <input type="checkbox"/> Oticon Real 1 miniRITE R | GTIN: (01) 05714464097302 |
| <input type="checkbox"/> Oticon Real 2 miniRITE R | GTIN: (01) 05714464097296 |
| <input type="checkbox"/> Oticon Real 3 miniRITE R | GTIN: (01) 05714464097289 |
| <input type="checkbox"/> Oticon Real 1 miniRITE T | GTIN: (01) 05714464097302 |
| <input type="checkbox"/> Oticon Real 2 miniRITE T | GTIN: (01) 05714464097296 |
| <input type="checkbox"/> Oticon Real 3 miniRITE T | GTIN: (01) 05714464097289 |

The following speakers are available for the above models:

- Speaker 60
- Speaker 85
- Speaker 100 (Power Instrument)
- Power flex mould, speaker 100 (Power Instrument)
- Power flex mould, speaker 105 (Power Instrument)
- MicroShell 85

## **Introduction to this booklet**

This booklet guides you on how to use and maintain your new hearing aids. Ensure you read this booklet carefully, including the Warnings section. This will help you get the most benefit from your hearing aids.

Your hearing care professional has adjusted the hearing aids to meet your needs. If you have additional questions, contact your hearing care professional.

A hearing care professional\* (hearing aid professional, audiologist, ENT (ear, nose and throat) doctor and hearing aid dispenser) is a person who is appropriately educated and has proven competency in professionally assessing hearing, selecting, fitting and delivering hearing instruments and rehabilitation care to persons with hearing loss. The education of the hearing care professional is in accordance with national or regional regulations.

\*The job title may vary from country to country.



## Intended use

<b>Intended use</b>	The hearing aid is intended to amplify and transmit sound to the ear.
<b>Indications for use</b>	Bilateral or unilateral impaired hearing of sensorineural, conductive or mixed type ranging from a slight (16 dB HL*) to profound (95 dB HL*) degree of hearing loss, with an individual frequency configuration.
<b>Intended user</b>	Person with hearing loss using a hearing aid and their caregivers. Hearing care professional responsible for adjusting the hearing aid.
<b>Intended user group</b>	Adults and children more than three years old..
<b>Use environment</b>	Indoor and outdoor.
<b>Contraindications</b>	Not suitable for infants below three years old. Users of active implants must pay special attention when using the hearing aid. For more information read the <b>Warnings</b> section.
<b>Clinical benefits</b>	The hearing aid is designed to provide better speech understanding to help ease communication with the aim of improving quality of life.

*\*As specified by the American Speech-Language-Hearing Association, [asha.org](http://asha.org), using pure-tone average of 0.5, 1 and 2 kHz.*

**IMPORTANT NOTICE**

The hearing aid amplification is uniquely adjusted and optimised to your personal hearing capabilities during the hearing aid fitting performed by your hearing care professional.

# Table of contents

## About

Your hearing aid, speaker and earpiece	11
--	----

## Getting started

Charging time	14
Battery performance	15
Identify left and right hearing aids	16
Turn hearing aids ON/OFF	18
Low battery indication	22
Replace disposable battery (size 312)	24

## Daily use

Put on the hearing aid	26
Cleaning	28
Clean the hearing aid	29
Replace the standard earpiece	30
ProWax miniFit filter	31
Replace ProWax miniFit filter	32
Clean customised earpiece	33
Replace ProWax filter	34



	Change volume	35
	Change program	36
	Store your hearing aids	38
<b>Options</b>		
	Flight mode	40
	Mute/unmute your hearing aids	42
	Use your hearing aids with an iPhone, iPad or iPod	44
	Use your hearing aids with Android devices	45
	Pairing and compatibility	46
	Call handling	47
	Wireless accessories and other options	48
<b>Tinnitus</b>		
	Tinnitus SoundSupport™ (optional)	50
	Guidelines for Tinnitus SoundSupport users	51
	Sound options and volume adjustments	52
	Limitation on use time	54
	Warnings related to tinnitus	57

Continues on next page

# Table of contents

## Warnings

---

General warnings	58
------------------	----

---

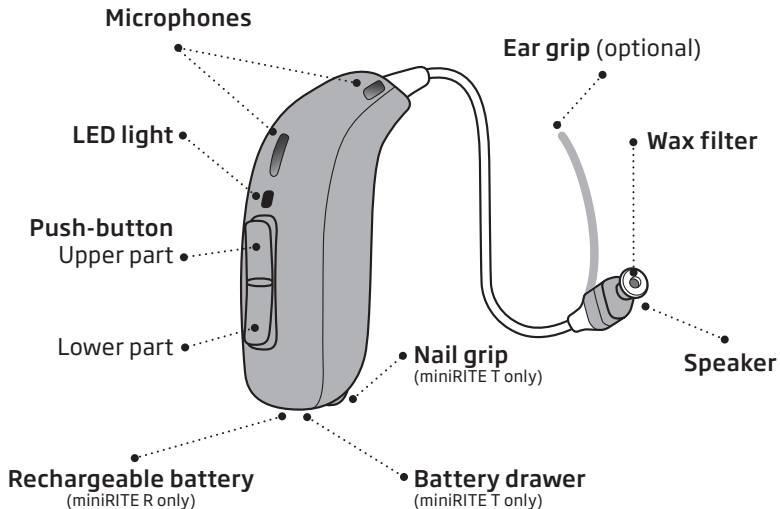
## More info

---

Troubleshooting	66
Water and dust resistant (IP68)	70
Conditions of use	71
Technical information	73
Australia only	80
Warranty	81
Your individual hearing aid settings	82
Sound and LED light indicators	85


---

# Your hearing aid, speaker and earpiece




The hearing aids use one of the following speakers:

### Standard speakers

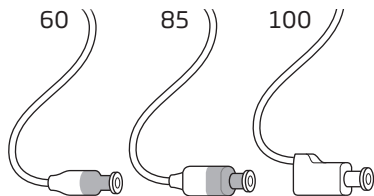

### Customised speakers

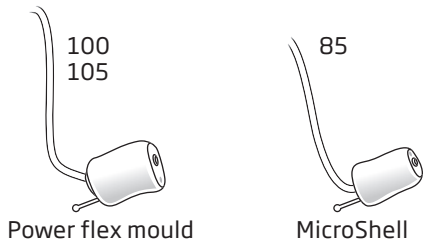

The speakers use one of the following earpieces:

### Standard earpieces



OpenBass dome



Bass dome, double vent



Power dome



Grip Tip

Available in small and large, left and right, with or without vent.

### Customised earpieces



LiteTip



Micro mould

### Dome sizes



5 mm\*



6 mm



8 mm



10 mm



12 mm

\*Only as OpenBass dome for speaker 60

### Note

For details on replacing the dome, see the **Replace standard earpieces** section.

## Charging time

(miniRITE R only)

Ensure you fully charge your hearing aids before first use and charge them every night. That ensures you start your day with fully charged hearing aids.

If your hearing aid's battery is completely drained, the normal charging time is:

3 hours	1 hour	0.5 hour
Fully charged	50% charged	25% charged

When the battery is fully charged, the charging process stops automatically.

Charging time may vary depending on the remaining capacity of the battery and between the left/right hearing aid.

For instructions on how to use your charger, see the charger's instructions for use.

## **Battery performance**

The battery performance varies depending on your individual use and hearing aid settings. Streaming sound from a TV, mobile phone or connectivity devices can influence this performance.

### **Rechargeable battery – miniRITE R only**

If your rechargeable hearing aids do not perform for a full day, contact your hearing care professional.

If your hearing aids run out of battery, ensure you recharge them by placing them in the charger.

Be aware that by restarting the hearing aids, it does not give you more usage time.

## Identify left and right hearing aids

It is important to distinguish between the left and the right hearing aid, as they may be programmed differently.

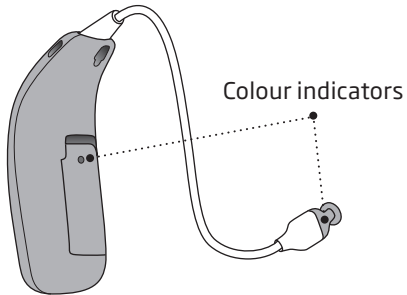
Colour indicators mark the left and right hearing aids. A **RED** indicator marks the RIGHT hearing aid. A **BLUE** indicator marks the LEFT hearing aid.

You can find colour indicators on the hearing aid itself and on 60 and 85 speakers as shown. They can also be found on 100 speakers and customised earpieces.


For 105 speakers, the colour indicator is found on the earpiece.



## miniRITE R



## miniRITE T

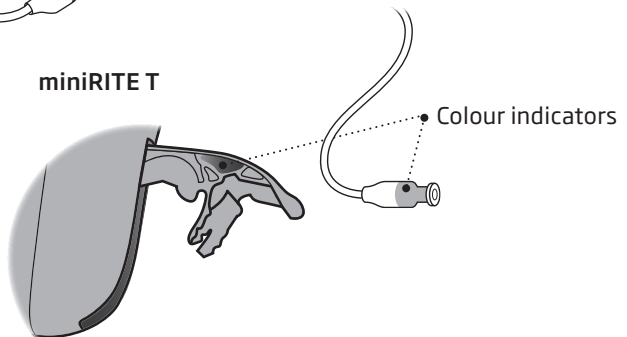

## Turn hearing aids ON/OFF

### Using the charger – miniRITE R only

Your hearing aids automatically turn ON when removed from the charger.\*

The hearing aid LED light turns **GREEN** after a few seconds. Wait until the hearing aid LED light blinks **GREEN** twice, confirming that it is ready for use. Depending on your hearing aid settings, you may also hear a start-up jingle.

Your hearing aid automatically turns OFF and starts charging when placed in the charger. The hearing aid LED light turns **ORANGE**.

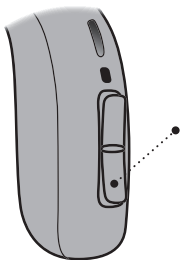
\* If the charger is disconnected from power or runs out of battery whilst charging your hearing aids, your hearing aids may turn OFF to save power. If your hearing aids do not automatically turn ON when removed from the charger, turn them on manually using the push-button.

**IMPORTANT NOTICE**

If applicable, ensure that your charger is powered or that the charger's built-in battery is charged when the hearing aid is seated in the charging port. For more information, see your charger's instructions for use.

## Using the push-button – miniRITE R only

Your hearing aids can be turned ON/OFF by using the push-button.



### To turn ON

Press and hold the lower part of the push-button for approximately two seconds until the hearing aid LED light turns **GREEN**.

Release the push-button and wait until the hearing aid LED light blinks **GREEN** twice.

The hearing aid is now turned ON.

### To turn OFF

Press and hold the lower part of the push-button for approximately three seconds until the hearing aid LED light turns **ORANGE**. The hearing aid plays four descending tones. Release the push-button and the hearing aid is turned OFF.

For information regarding tones, see the **Sound and LED light indicators** section.

## Using the battery drawer – miniRITE T only

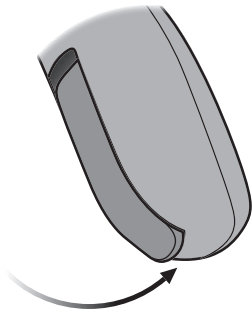
The battery drawer is used to turn the hearing aids ON and OFF. Depending on your hearing aid settings, you may also hear a start-up jingle.

### LED light

Three green blinks indicate that the hearing aid is turned on.


#### Turn ON

Close the battery drawer with the battery in place.



#### Turn OFF

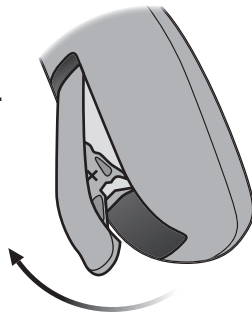Open the battery drawer.



To save battery life, make sure your hearing aids are switched OFF when you are not wearing them. To perform a quick reset of hearing aid settings, open and close the battery drawer.


## Low battery indication

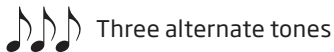Your hearing aids indicate when the battery is running low either by playing three alternate tones or a Voice Notification (optional) in selected languages. The indication is repeated periodically until the battery runs out. Your hearing care professional can set your hearing aids to match your preferences. To extend battery performance, ensure you stop any audio streaming.

- miniRITE R: Your hearing aid will indicate that the battery is running low. This gives you approximately two hours before the hearing aid runs out of battery. At this point, you may continue to stream audio for approximately one hour.
- miniRITE T: Your hearing aid will indicate that the battery is running low. This gives you approximately 15 minutes before the hearing aid runs out of battery. At this point, Bluetooth® connectivity is turned OFF.


Just before the battery runs out completely, you hear four descending tones.

## The battery is running low



Three alternate tones

## The battery has run out

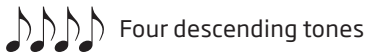

Four descending tones



Voice Notification


## Hearing aid LED light

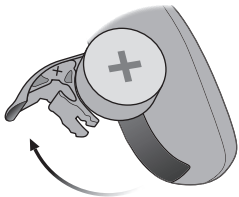Continuous **ORANGE** blinks indicate low battery.

# Replace disposable battery (size 312)


(miniRITE T only)

## 1. Remove



Fully open the battery drawer and remove the battery.

## 2. Uncover

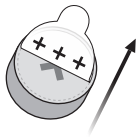


Remove the sticky label from the + side of the new battery.

### Tip:

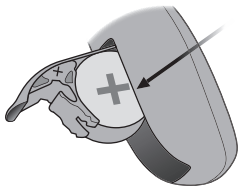Wait two minutes so that the battery can draw air, to ensure optimal functioning.

## 3. Insert




Insert the new battery into the battery drawer. Insert it from above, NOT from the side.

Ensure that the + side of the battery and battery drawer align.




#### 4. Close

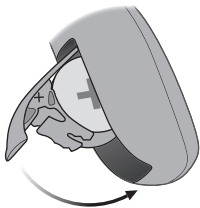

Close the battery drawer.  
You may hear a jingle through the  
earpiece.

#### Tip

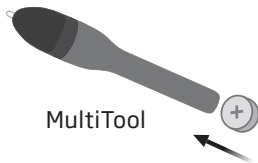

You can use the MultiTool to change  
the battery. Use the magnetic end to  
remove and insert batteries.

# Put on the hearing aid

## Step 1

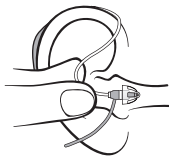


Place the hearing aid behind your ear.

You should always use the speaker with an earpiece attached.

Use only parts designed for your hearing aid.


## Step 2



Hold the bend of the speaker wire between your thumb and index finger.

The earpiece should point towards the opening of the ear canal.

## Step 3

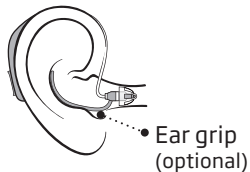


Gently push the earpiece into your ear canal until the speaker wire sits close to your head.

If the speaker has an ear grip, place it in your ear so it follows the contour of the ear.



## Cleaning

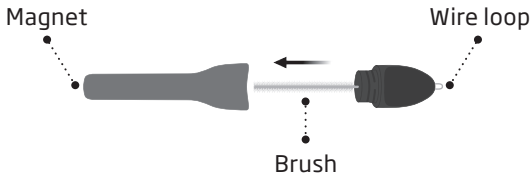The MultiTool contains a brush and wire loop for cleaning and removing earwax. If you need a new MultiTool, contact your hearing care professional.



### IMPORTANT NOTICE

The MultiTool has a built-in magnet. Keep the MultiTool at least 30 centimetres (one foot) away from credit cards and other magnetically sensitive devices.

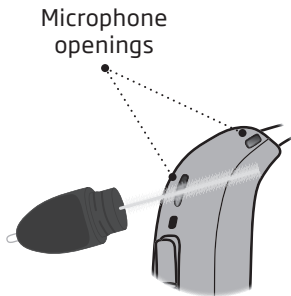
## Clean the hearing aid

When handling the hearing aid, hold it over a soft surface to avoid damage in case you drop it.

### Clean the microphone openings

Use the brush of the MultiTool to carefully brush debris away from the openings and the surface around the openings.

Ensure that you do not forcefully squeeze parts of the MultiTool into the microphone openings. This may damage the hearing aid.




### IMPORTANT NOTICE

To clean the hearing aids, use a soft, dry cloth. The hearing aids must never be washed or immersed in water or other liquids.

## Replace the standard earpiece


The standard earpiece (dome or Grip Tip) should not be cleaned. If the earpiece is filled with earwax, replace it with a new one. Grip Tips should be replaced at least once a month.
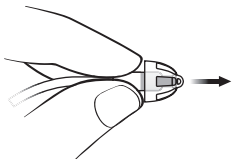
### Step 1




Hold on to the speaker and pull off the earpiece.

### Step 2

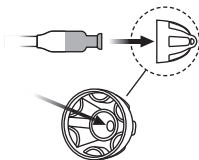

Insert the speaker exactly into the middle of the earpiece to get a secure fit.

### Step 3

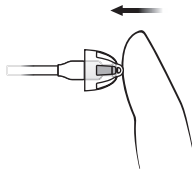


Push firmly to ensure that the earpiece is securely fastened.

## ProWax miniFit filter

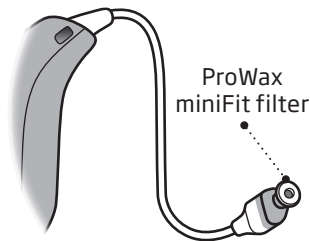The speaker has a white wax filter attached to the end where the earpiece is attached. The wax filter keeps earwax and debris from damaging the speaker.

Ensure you replace the filter when clogged, or if the hearing aid does not sound normal. Alternatively, contact your hearing care professional. Ensure you remove the earpiece from the speaker before replacing the wax filter. To do this, see the **Replace standard earpieces** section.




### IMPORTANT NOTICE

Ensure you always use the same type of wax filter as was originally supplied with the hearing aids. If you are in doubt about the use or replacement of wax filters, contact your hearing care professional.


# Replace ProWax miniFit filter
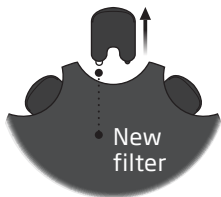
## 1. Tool




Remove the tool from the shell. The tool has two pins: one empty for removal and one with the new wax filter.

## 2. Remove

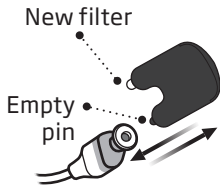

Insert the empty pin into the wax filter in the speaker and pull it out.

## 3. Insert

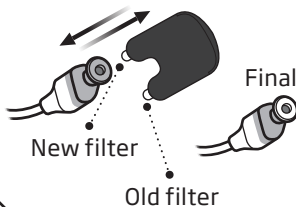

Insert the new wax filter using the other pin. Remove the tool and throw it out.

### Note

If you use a mould or LiteTip, your hearing care professional must replace the wax filter in the speaker.



## Clean customised earpiece

The earpiece should be cleaned regularly.

The earpiece has a white wax filter that keeps earwax and debris from damaging the speaker. Ensure you replace the filter when clogged, or if the hearing aid does not sound normal.

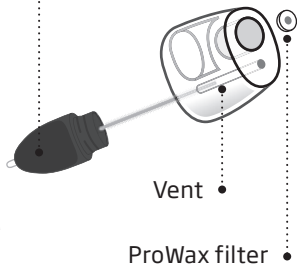
Alternatively, contact your hearing care professional.

### Note

If you use a mould or LiteTip, your hearing care professional must replace the wax filter in the speaker.

\* LiteTip does not have a wax filter.

- Clean the vent by inserting the brush through the hole, twisting it slightly.




### IMPORTANT NOTICE

If the earpiece is not on the speaker when removed from the ear, the earpiece may still be in the ear canal. For further instructions, consult your hearing care professional.


# Replace ProWax filter
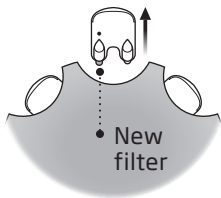
## 1. Tool




Remove the tool from the shell. The tool has two pins: one empty for removal and one with the new wax filter.

## 2. Remove

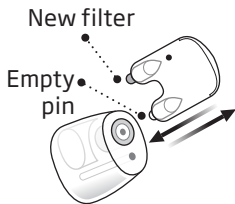

Insert the empty pin into the wax filter in the earpiece and pull it out.

## 3. Insert

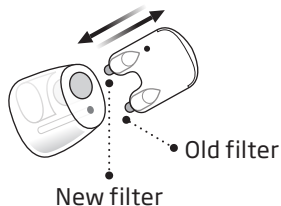


Insert the new wax filter using the other pin. Remove the tool and throw it out.

## Change volume

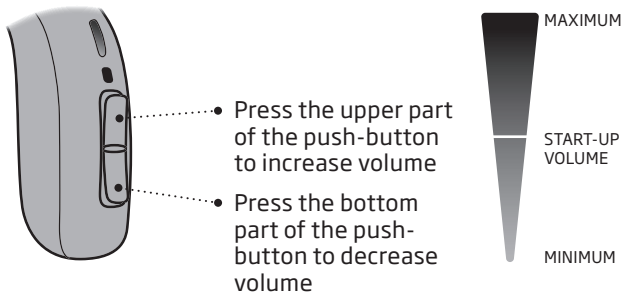
The push-button allows you to adjust the volume. When you increase or decrease the volume, you hear a beep.

For information regarding button press times, see the table **General settings overview for your hearing aid, in the Your individual hearing aid settings** section at the end of this booklet.



## Change program

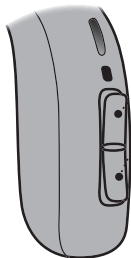Your hearing aid can have up to four different programs, which are programmed by your hearing care professional. Depending on the program, you will hear one to four tones when you change the program.



- Press the push-button to switch between programs.

The program cycle switches one program forward when the upper part of the push-button is pressed, for example program one to two or program four to one.

If the lower part of the push-button is pressed, the program cycle goes backward, for example two to one or program one to four.

For information regarding tones, see the **Sound and LED light indicators** section.

For information regarding button press times, see the table **General settings overview for your hearing aid**, in the **Your individual hearing aid settings** section at the end of this booklet.

# Store your hearing aids

(miniRITE R only)

When you are not using your hearing aids, the charger is the best place to keep them.

To ensure the longest life of the rechargeable battery in your hearing aids, do not expose them to excessive heat. For example, do not leave your hearing aids in the sun in front of a window or in a car, even if they are in the charger.

## Long-term storage

Before you put away or store your hearing aids for a prolonged period (more than 14 days), ensure you first fully charge them, and then turn them OFF. This way the battery can be charged again.

### Note

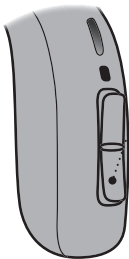
To protect the rechargeable battery, it is necessary that you fully charge the hearing aid every six months. If a stored hearing aid is not charged within a six-month period, the rechargeable battery must be replaced.



## Flight mode

When Flight mode is activated, Bluetooth connectivity is turned OFF. However, your hearing aids are still turned ON and functioning. Be aware that pressing the push-button on one hearing aid, activates Flight mode on both hearing aids. For information about sounds and lights, see the **Sound and LED light indicators** section.

miniRITE R



- **To activate and deactivate**

Press and hold the lower part of the push-button for approximately seven seconds.

Four descending tones, a jingle and an LED light pattern will confirm your action.




**To activate and deactivate**

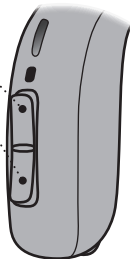Press and hold either part of the push-button for approximately seven seconds.

A jingle and an LED light pattern will confirm your action.

Opening and closing the battery drawer also deactivates Flight mode.

**miniRITE T**



## Mute/unmute your hearing aids

You can mute your hearing aids by using one of the following optional app and devices:

- Oticon Companion
- ConnectClip
- Remote Control 3.0

### How to unmute your hearing aids

You can unmute your hearing aids by using one of the optional app and devices or by applying a short press to the upper or lower part of the push-button on the hearing aids.

### miniRITE T only

You can also mute your hearing aids by pressing either part of the push-button for approximately four seconds.

**IMPORTANT NOTICE**

Do not use the mute function as an OFF button, as the hearing aid is still using battery power in this mode.

## Use your hearing aids with an iPhone, iPad or iPod

Your hearing aids are Made for iPhone hearing aids and allow for direct streaming from your iPhone, iPad or iPod touch.\*

Oticon Companion can be used to control your hearing aids from your mobile device.\*

Additionally, your hearing aids offer hands-free communication with compatible iPhone and iPad devices.\*

For assistance in using your hearing aids with any of these products, contact your hearing care professional.

Apple, the Apple logo, iPhone, iPad and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

## **Use your hearing aids with Android devices**

Your hearing aids support Audio Streaming for Hearing Aids (ASHA) and allow for direct streaming from selected Android™ devices.\*

Oticon Companion can be used to control your hearing aids from your mobile device.\*

For assistance in using your hearing aid with Android devices, contact your hearing care professional.

\* For a list of compatible iPhone, iPad, iPod touch and Android devices, visit: [www.oticon.com.au/compatibility](http://www.oticon.com.au/compatibility).

## Pairing and compatibility

For instructions on how to pair your hearing aids with iPhone, iPad, iPod touch or Android devices, visit:

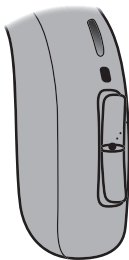
[www.oticon.com.au/pairing](http://www.oticon.com.au/pairing)

For a list of compatible iPhone, iPad, iPod touch and Android devices, visit:

[www.oticon.com.au/compatibility](http://www.oticon.com.au/compatibility)

## Call handling

You can answer, reject or end phone calls with the push-buttons on your hearing aids. For this functionality, your hearing aids must be paired with a compatible iPhone.\*



- **To accept**

Use a short press on either part of the push-button to accept a phone call. A short tone confirms your action.

**To reject**

Press and hold on either part of the push-button to reject a phone call. Short, descending tones will confirm your action.

**To end**

Press and hold on either part of the push-button to end a phone call. Short, descending tones confirm your action.

\* For a list of compatible iPhones, visit: [www.oticon.com.au/compatibility](http://www.oticon.com.au/compatibility).

## Wireless accessories and other options

There are a range of accessories available as an enhancement to your wireless hearing aids. These enable you to hear and communicate better in everyday situations.

### ConnectClip

A device that can be used as a remote microphone and hands-free headset when paired to your mobile phone.

### Phone Adapter 2.0

A device that when used together with hearing aids and ConnectClip, lets you communicate hands-free via a landline phone.

### TV Adapter 3.0

A device that streams sound from a TV or electronic audio device to your hearing aids.

### Remote Control 3.0

A device that lets you change program, adjust volume or mute your hearing aids.



## **EduMic**

A device that can be used as a remote microphone in classrooms, work environments, public places (using Telecoil) and other settings.

## **Oticon Companion**

An application that gives you wireless access to your hearing aids through a mobile device for control of functionality, and that allows an update of your hearing aids remotely. For iPhone, iPad, iPod touch and Android devices. Ensure that you only download Oticon Companion from the official app stores.

Apple, the Apple logo, iPhone, iPad and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, Google Play and the Google Play logo are trademarks of Google LLC.

## **Telecoil**

Telecoil can help you hear better when using a phone with a built-in loop or when in buildings with teleloop systems such as theatres, places of worship or lecture rooms. This symbol is shown wherever a teleloop has been installed.



## **Tinnitus SoundSupport™ (optional)**

### **Intended use of Tinnitus SoundSupport**

Tinnitus SoundSupport is a tool intended to generate sounds to provide temporary relief for patients suffering from tinnitus as part of a tinnitus management program.

The target population is the adult population over 18 years of age.

Tinnitus SoundSupport is targeted to licenced hearing care professionals (audiologists, hearing aid specialists, or otolaryngologists) who are familiar with the evaluation and treatment of tinnitus and hearing loss. Fitting of Tinnitus SoundSupport must be done by a hearing care professional participating in a tinnitus management program.

## **Guidelines for Tinnitus SoundSupport users**

These instructions contain information about Tinnitus SoundSupport, which may have been enabled in your hearing aids by your hearing care professional.

Tinnitus SoundSupport is a tinnitus management device intended to generate sound of sufficient intensity and bandwidth to help manage tinnitus.

Your hearing care professional will also be able to offer the appropriate follow-up care. It is important to follow his/her advice and directions regarding such care.

### **Prescription use only**

Good health practice requires that a person reporting tinnitus have a medical evaluation by a licenced ear doctor before using a sound generator. The purpose of such an evaluation is to ensure that any medically treatable condition that may cause tinnitus is identified and treated prior to using a sound generator.

## Sound options and volume adjustments

Tinnitus SoundSupport is programmed by your hearing care professional to match your hearing loss and preferences for tinnitus relief. It offers several different sound options. Together with your hearing care professional, you can select your preferred sound(s).

### Tinnitus SoundSupport programs

Together with your hearing care professional, you decide which programs you may want to have Tinnitus SoundSupport activated for. The sound generator can be activated in up to four different programs.

### Mute

If you are in a program for which Tinnitus SoundSupport is activated, the mute functionality only mutes the environmental sounds, and not the sound from Tinnitus SoundSupport. For information on how to mute your hearing aids, see the **Mute your hearing aids** section.

## **Volume adjustments with Tinnitus SoundSupport**

When you select a hearing aid program for which Tinnitus SoundSupport is activated, your hearing care professional can only set the push-button on your hearing aid to work as a volume control for the tinnitus relief sound.

Your hearing care professional sets the volume control for the sound generator in one of two ways:

- A) Change volume in each ear separately, or
- B) Change volume in both ears simultaneously.

For more information about volume adjustments with Tinnitus SoundSupport, see the table **Tinnitus SoundSupport settings overview for your hearing aid** in the **Your individual hearing aid settings** section at the end of this booklet.

## Limitation on use time

### Daily use

The volume levels of Tinnitus SoundSupport can be set to a level which could lead to permanent hearing damage when used for a prolonged period. Your hearing care professional will advise you of the maximum amount of time per day you should use Tinnitus SoundSupport. It should never be used at uncomfortable levels.

See the table **Tinnitus SoundSupport: Limitation on use**, in the **Your individual hearing aid settings** section at the end of this booklet to learn how many hours per day you can safely use the relief sound in your hearing aid.

# Important information for hearing care professionals about Tinnitus SoundSupport

## Device description

Tinnitus SoundSupport is a module function that can be enabled in the hearing aids by the hearing care professional.

## Maximum wearing time

The wearing time of Tinnitus SoundSupport will decrease as you increase the level above 80 dB(A) SPL. The fitting software automatically displays a warning when the hearing aid exceeds 80 dB(A) SPL. See **Max wearing time indicator** next to the tinnitus fitting graph in the fitting software.

## The volume control is deactivated

By default, the volume control for the sound generator is deactivated in the hearing aid. Risk of noise exposure increases when the volume control is activated.

## If the volume control is activated

A warning may be displayed if you activate the tinnitus volume control in the **Buttons & Indicators** screen. This occurs if the relief sound can be listened to at levels that may cause hearing damage.

The **Max wearing time** table in the fitting software displays the number of hours the patient can safely use Tinnitus SoundSupport.

1. Note the max wearing time for each program for which Tinnitus SoundSupport is activated.
2. Write those values in the table **Tinnitus SoundSupport: Limitation on use**, in the back of this booklet.
3. Instruct your patient accordingly.



## **Warnings related to tinnitus**

If your hearing care professional has activated the sound generator Tinnitus SoundSupport, please pay attention to the following warnings.

There are some potential concerns associated with the use of any sound generated by a tinnitus management device. Among them are the potential worsening of tinnitus, and/or a possible change in hearing thresholds.

Should you experience or notice a change in hearing or tinnitus, or any dizziness, nausea, headaches, heart palpitations or possible skin irritation at the point of contact with the device, you should immediately discontinue use of the device and consult a medical, audiology or other hearing care professional.

As with any device, misuse of the sound generator feature may cause potentially harmful effects. Care should be taken to prevent unauthorised use and to keep the device out of reach of children and pets.

### **Maximum wearing time**

Always follow the maximum wearing time per day of the Tinnitus SoundSupport advised by your hearing care professional. Prolonged use may lead to worsening of your tinnitus or of your hearing loss.

## **General warnings**

For your personal safety and to ensure correct usage, you should familiarise yourself fully with the following general warnings before using your hearing aids.

Consult your hearing care professional if you experience unexpected operations or serious incidents with your hearing aids during use or because of its use. Your hearing care professional will support you with issue handling and, if relevant, reporting to the manufacturer and/or the national authorities.

Note that hearing aids do not restore normal hearing and do not prevent or improve a hearing impairment resulting from organic conditions. Hearing aids are only a part of hearing rehabilitation and may need to be supplemented by auditory training and instruction in lipreading. Furthermore in most cases, infrequent use of hearing aids results in a user not attaining full benefit from it.

(miniRITE R only)

Only charge the hearing aids with a designated charger. Other chargers risk destroying the hearing aids and batteries.

This hearing aid is supported by a non-removable rechargeable lithium-ion battery cell. Please ensure to charge the hearing aid and familiarise yourself with the safety and handling information related to rechargeable hearing aids.

Do not try to get access to the battery inserted in the hearing instrument. The battery must only be replaced by your hearing care professional.

### **Usage of hearing aids**

Hearing aids should be used only as directed and adjusted by your hearing care professional. Misuse can result in sudden and permanent hearing loss.

Never allow others to wear your hearing aid, as incorrect usage could cause permanent damage to their hearing.

### **Choking hazards and risk of swallowing batteries or other small parts**

Hearing aids, their parts and batteries should be kept out of reach of children and anyone who might swallow these items or otherwise cause injury to themselves.

*If a battery, hearing aid or small part is swallowed, see a doctor immediately.*

(miniRITE T only)

Batteries have occasionally been mistaken for pills. Therefore, check your medicine carefully before swallowing any pills.

### **Battery use (miniRITE T only)**

Always use batteries recommended by your hearing care professional. Batteries of low quality may leak and cause bodily harm.

Never attempt to recharge your batteries and never dispose of batteries by burning them. There is a risk that the batteries will explode.

### **Explosives (miniRITE R only)**

Your hearing aids are safe to use under normal usage conditions. The hearing aid has not been tested for compliance with international standards concerning explosive environments.

Therefore, do not use the hearing aid in environments with danger of explosions (e.g., mines, oxygen rich environments or areas where flammable anaesthetics are handled).

### **Fatality hazards and risk of swallowing lithium-ion batteries or placing them in the ear or nose (miniRITE R only)**

Never swallow lithium-ion batteries or place them in the ear or nose, as this may lead to serious injury or death in as little as two hours. This can be due to chemical burns, which can permanently damage the nose or ear, or potentially lead to perforation of the inner organs.

If a lithium-ion battery is swallowed or

Continues on next page

## **General warnings**

placed in the ear or nose, seek emergency medical treatment immediately. Keep the batteries in the original packaging until use. Dispose of used batteries immediately.

### **Rechargeable battery** (miniRITE R only)

Do not attempt to open the hearing aid, as it may damage the battery.

Never attempt to replace the battery. If battery replacement is needed, please return your device to the supplier. The service guarantee is void if there are signs of tampering.

In case of battery leakage do not wear or carry your hearing aid, as it may cause skin irritation due to acids leaking from the battery. If your skin has been in contact with the leaked battery acids, use a wet cloth to wipe it off and ensure no acid is left on your skin. If you experience skin irritation, consult your doctor. For further handling instructions of your hearing aid consult your hearing care professional.

The safety of recharging batteries using a USB connector is determined by the external signal source. When connected to external equipment plugged into a power point, this equipment must comply with IEC 62368-1 or equivalent safety standards.

### **Dysfunction**

Be aware of the possibility that your hearing aids may stop working without notice. Keep this in mind when you depend on warning sounds (e.g., when you are in traffic). The hearing aids may stop functioning, for instance if the batteries have expired or if the tubing is blocked by moisture or earwax.

## **Active implants**

The hearing aids have been thoroughly tested and characterised for human health according to international standards for human exposure (Specific Absorption Ratio – SAR), induced electromagnetic power and voltages into the human body.

The exposure values are well below internationally accepted safety limits for SAR, induced electromagnetic power and voltages into the human body defined in the standards for human health and coexistence with active medical implants such as pacemakers and heart defibrillators.

If you have an active brain implant, please contact the manufacturer of your implantable device for information about the risk of disturbance.

The Autophone magnet or MultiTool (which has a built-in magnet) should be kept more than 30 centimetres (1 foot) away from the implant, e.g., do not carry it in your breast pocket.

Follow the guidelines recommended by the manufacturers of implantable defibrillators and pacemakers regarding their use with magnets.

## **Cochlear implants**

If you are using a cochlear implant (CI) on one ear and a hearing aid on the other ear, make sure to always keep your CI at least one centimetre (0.4 inches) from your hearing aid. The magnetic field from CI sound processors, coils and magnets may permanently damage the speaker unit in your hearing aid. Never place the devices close together on a table (e.g., when cleaning or changing batteries). Do not carry the CI and the hearing aid together in the same box.

Continues on next page

## **General warnings**

### **Detached earpiece in the ear canal**

If the earpiece is not on the speaker when removed from the ear, the earpiece may still be in the ear canal. For further instructions, consult your hearing care professional.

### **X-ray/CT/MR/PET scanning, electrotherapy and surgery**

Remove your hearing aid before X-ray examinations and CT/MR/PET scans, electrotherapy, surgery, etc. as your hearing aid may be damaged when exposed to strong electromagnetic fields.

### **Heat and chemicals**

Your hearing aids must never be exposed to extreme heat (e.g., left inside a parked car in the sun).

The hearing aids must not be dried in microwave ovens or other ovens.

The chemicals in cosmetics, hairspray, perfume, aftershave lotion, sunscreen lotion and insect repellent can damage the hearing aids. Always remove your hearing aids before applying such products and allow time to dry before use.

### **Connection to external equipment**

The safety of the hearing aids, when connected to external equipment with USB cable and/or directly, is determined by the external signal source. When the hearing aids are connected to external equipment which is plugged into a power point, this equipment must comply with IEC 62368-1 or equivalent safety standards.

### **Power hearing aid**

Special care should be exercised in selecting, fitting, and using hearing aids where the maximum sound pressure capability exceeds 132 dB SPL as there may be a risk of impairing the remaining hearing of the hearing aids user.

*For information on whether your instrument is a power hearing aid, see the model overview section in this booklet.*

### **Possible side effects**

Hearing aids, moulds or domes may cause an accelerated accumulation of earwax.

The non-allergenic materials used in hearing aids may, in rare cases cause a skin irritation or other side effects.

*If these conditions occur, seek consultation with a doctor.*

### **Use on aircraft**

Your hearing aids have Bluetooth wireless technology. On board an aircraft, the hearing aids must be put into Flight mode to deactivate Bluetooth, unless Bluetooth is permitted by the flight personnel.

### **Use of third-party accessories**

Only use accessories, transducers or cables supplied by the manufacturer. Non-original accessories may result in reduced electromagnetic compatibility (EMC) of your hearing aids.

### **Modification of hearing aids is not allowed**

Changes or modifications not expressly approved by the manufacturer will void the warranty of the equipment.

Continues on next page

## **General warnings**

### **Interference**

The hearing aids have been thoroughly tested for interference according to the most stringent international standards, including EN/IEC 60601-1-2 and its amendments. Electromagnetic interference may occur in the vicinity of equipment with the symbol to the left. Portable and mobile RF (radio frequency) communications equipment can affect the performance of your hearing aids. If your hearing aids are affected by electromagnetic interference, move away from the source to reduce the interference.



### **WARNING**

New and used button/coin batteries are hazardous and should be kept away from children. Button/coin batteries can cause severe or fatal injuries in 2 hours or less if swallowed or placed inside any part of the body. Seek medical attention immediately if this is suspected.

If it is suspected a button/coin battery has been swallowed or otherwise placed inside any part of the body, you should contact the

Australian Poisons Information Centre  
on 13 11 26

New Zealand National Poisons Centre  
on 0800 764 766  
immediately for 24/7 fast, expert advice.





# Troubleshooting

## miniRITE R and miniRITE T

Symptom	Possible causes
<b>No sound</b>	Hearing aid is out of power
	Dead battery
	Clogged earpieces (dome, Grip Tip, or mould)
	Hearing aid microphone muted
<b>Intermittent or reduced sound</b>	Clogged sound outlet
	Moisture
<b>Squealing noise</b>	Hearing aid earpiece incorrectly inserted
	Earwax accumulated in ear canal
<b>Beeping</b>	If your hearing aid plays eight beeps, four times consecutively, your hearing aid needs a microphone service check
<b>Pairing issue with smartphone</b>	Bluetooth connection failed
	Only one hearing aid is paired

## Solutions

Charge the hearing aid (**miniRITE R only**) / Replace the battery (**miniRITE T only**)

Contact your hearing care professional (**miniRITE R only**) / Replace the battery (**miniRITE T only**)

Clean mould

Replace wax filter, dome or Grip Tip

Unmute the hearing aid microphone

Clean mould or replace wax filter, dome or Grip Tip

Gently wipe the hearing aid and let it dry

Re-insert the earpiece

Have ear canal examined by your doctor

Contact your hearing care professional

1) Unpair your hearing aid

2) On your phone, turn Bluetooth OFF and ON again

3) Turn the hearing aid OFF and then turn it back ON

4) Pair your hearing aid again (for guidance, visit: [www.oticon.com.au/pairing](http://www.oticon.com.au/pairing))

### Note

If none of the above solutions work, consult your hearing care professional.

Continues on next page

# Troubleshooting

**miniRITE R only** – To troubleshoot the charger, see your charger's instructions for use.

Symptom	Possible causes
<b>The hearing aid LED light remains turned OFF when the hearing aid is placed in the charger</b>	The charger is not turned ON
	The hearing aid or charger's temperature is either too warm or too cold
	Charging is incomplete. The charger has stopped charging to protect the battery.
	The hearing aid is incorrectly seated in the charger
<b>The hearing aid LED light blinks GREEN when the hearing aid is placed in the charger</b>	The hearing aid has not been in use for a prolonged period
<b>The hearing aid LED light blinks ORANGE when the hearing aid is placed in the charger</b>	System error

## Solutions

Verify that the charger's power plug is correctly connected, or the power bank has enough battery

Move the charger and hearing aid to a location with a temperature between +5°C and +40°C (+41°F and +104°F)

Reinsert the hearing aid into the charger. This completes the charging within approximately 15 minutes.

Check the charging ports for foreign objects

Depending on how depleted the battery in the hearing aid is, the hearing aid automatically resumes normal charging mode after a timeframe of up to 10 minutes. Ensure you leave the hearing aids in the charger during this process.

Contact your hearing care professional

### Note

If none of the above solutions work, consult your hearing care professional.

## Water and dust resistant (IP68)

Your hearing aid is dust tight and protected against the ingress of water, which means it is designed to be worn in all daily life situations. The water and dust resistance means you do not have to worry about your hearing aid getting wet when it rains or comes into contact with sweat.

### **miniRITE R only**

Before charging the hearing aid, make sure to wipe off any moisture.

### **miniRITE T only**

Should your hearing aid come into contact with water and stop working, follow these guidelines:

1. Gently wipe off any water
2. Open the battery drawer, remove the battery and gently wipe off any water in the battery drawer
3. Let the hearing aid dry with the battery drawer left open for approximately 30 minutes
4. Insert a new battery.

### **IMPORTANT NOTICE**

Do not wear your hearing aid while showering or participating in water activities. Do not immerse your hearing aid in water or other liquids.

# Conditions of use

(miniRITE R only)

<b>Operating conditions</b>	Temperature: +5°C to +40°C (41°F to 104°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa
<b>Charging conditions</b>	Temperature: +5°C to +40°C (41°F to 104°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa
<b>Transportation and storage conditions</b>	Temperature and humidity should not exceed the below limits for extended periods during transportation and storage:  <b>Transportation:</b> Temperature: -20°C to +60°C (-4°F to 140°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa <b>Storage:</b> Temperature: -20°C to +30°C (-4°F to 86°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa

## Note

For more information about the charger's conditions of use, see your charger's instructions for use.

## Conditions of use

(miniRITE T only)

<b>Operating conditions</b>	Temperature: +1°C to +40°C (34°F to 104°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa
<b>Transportation and storage conditions</b>	Temperature and humidity should not exceed the below limits for extended periods during transportation and storage:  <b>Transportation:</b> Temperature: -25°C to +60°C (-13°F to 140°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa <b>Storage:</b> Temperature: -25°C to +60°C (-13°F to 140°F) Humidity: 5% to 93% relative humidity, non-condensing Atmospheric pressure: 700 hPa to 1060 hPa



## Technical information

The hearing aid contains the following two radio technologies:

The hearing aid contains a radio transceiver using short range magnetic induction technology operating at 3.84 MHz. The magnetic field strength of the transmitter is very weak and always below 15 nW (typically below  $-40 \text{ dB}\mu\text{A}/\text{m}$  ( $-12.20 \text{ dB}\mu\text{A}/\text{ft}$ ) at 10 metres (33 feet) distance).

The hearing aid also contains a radio transceiver using Bluetooth Low Energy technology and a proprietary short-range radio technology, both operating at ISM band 2.4 GHz.

The radio transmitter is weak and always below 9 mW equal to 9.6 dBm in total radiated power.

Only use your hearing aids in areas where wireless transmission is permitted.

The hearing aid complies with international standards concerning radio transmitters, electromagnetic compatibility and human exposure.

Due to the limited space available on the hearing aid, relevant approval markings can be found in this booklet. Additional information can be found in the Technical Data Sheet on [www.oticon.com.au](http://www.oticon.com.au)

The manufacturer declares that this hearing aid is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

This medical device complies with Medical Device Regulation (EU) 2017/745.

Declaration of Conformity is available from the headquarters.

Oticon A/S  
Kongebakken 9  
DK-2765 Smørum  
Denmark  
[www.oticon.global/doc](http://www.oticon.global/doc)

CE 0123



SBO Hearing A/S  
Kongebakken 9  
DK-2765 Smørum  
Denmark



Waste from electronic equipment must be handled according to local regulations.



 **Bluetooth®**

**IP68**

**Australian Sponsor**

Audmet Australia Pty Ltd  
Suite 4, Level 4, Building B  
11 Talavera Road  
North Ryde NSW 2113

## Description of symbols and abbreviations used in this booklet



### Warnings

Text marked with a warning symbol must be read before using the device.



### Manufacturer

The device is produced by the manufacturer whose name and address are stated next to the symbol. Indicates the medical device manufacturer, as defined in EU Regulations 2017/745.



### CE mark

The device complies with all required EU regulations and directives. The four-digit number indicates the identification of the notified body.



### Electronic waste (WEEE)

Recycle hearing aids, accessories or batteries according to local regulations. Hearing aid users can also return electronic waste to their hearing care professional for disposal. Electronic equipment covered by Directive 2012/19/EU on waste and electrical equipment (WEEE).



### Regulatory Compliance Mark (RCM)

The device complies with electrical safety, EMC and radio spectrum requirements for devices supplied to the Australian or New Zealand markets.

IP68

### IP code

Indicates the class of protections against harmful ingress of water and particulate matter according to EN 60529. IP6X indicates total dust protection. IPX8 indicates the protection against the effects of continuous immersion in water.



### Bluetooth logo

Registered trademark of Bluetooth SIG, Inc. where any use of such requires a licence.

## Description of symbols and abbreviations used in this booklet



### **Made for Apple badges**

Indicates that the device is compatible with iPhone, iPad and iPod touch.



### **Android badge**

Indicates that the device is compatible with Android.



### **Hearing loop**

This logo incorporates the universal symbol for hearing assistance. The "T" signifies that a hearing loop is installed.



### **Radio Frequency (RF) transmitter**

Your hearing aid contains an RF transmitter.

**GTIN**

### **Global Trade Item Number**

A globally unique 14-digit number used to identify medical device products including medical device software. GTIN in this booklet is related to medical device firmware (FW). GTIN on the regulatory packaging label is related to medical device hardware.

**FW**

### **FW**

Firmware version used in the device.



### **Button Battery Warning**

The device is intended to be used with button batteries. Refer to enclosed warnings before using the device.

## Description of additional symbols used on labels



### **Keep dry**

Indicates a medical device that needs to be protected from moisture.



### **Catalogue number**

Indicates the manufacturer's catalogue number so that the medical device can be identified.



### **Serial number**

Indicates the manufacturer's serial number so that a specific medical device can be identified.



### **Medical Device**

The device is a medical device.



Li-ion

### **Battery recycling symbol**

Li-ion battery recycling symbol.



### **Temperature limit**

Indicates the temperature limits to which the medical device can be safely exposed.



### **Humidity limitation**

Indicates the range of humidity to which the medical device can be safely exposed.



### **Radio Frequency Identification**

Indicates the presence of a passive radio-frequency identification tag incorporated into the device for manufacturing and service purposes.



### **Unique device identifier**

Indicates a carrier that contains unique device identifier information

## International warranty

Your device is covered by an international limited warranty issued by the manufacturer for a period of 12 months from the date of delivery to you. This limited warranty covers manufacturing and material defects in the device itself, but not in accessories such as batteries, tubing, speakers, earpieces and filters, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, physical changes in your ear, damage due to foreign objects entering the device, or incorrect adjustments are NOT covered by the limited warranty and may void it. Under this warranty, the manufacturer will choose, at its sole discretion, whether to repair the device, or replace it with

an equivalent model. The above warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods. Your hearing care professional may have issued a warranty that goes beyond the clauses of this limited warranty. Please consult him/her for further information.

### **If you need service**

Take your device to your hearing care professional, who may be able to sort out minor problems and make adjustments immediately. Your hearing care professional can assist you to obtain warranty service from the manufacturer. Your hearing care professional may charge a fee for their services.

## Australia only

The rights given under this warranty are in addition to other rights and remedies under a law in relation to the good or services to which the warranty relates to. Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Australian customers should contact Oticon Australia at Suite 4, Level 4, Building B, 11 Talavera Road, North Ryde NSW 2113. Telephone: 1800 633 738 if you have any questions regarding this warranty.



# Warranty

## Certificate

Name of owner: \_\_\_\_\_

Hearing care professional: \_\_\_\_\_

Hearing care professional's address: \_\_\_\_\_

Hearing care professional's phone: \_\_\_\_\_

Purchase date: \_\_\_\_\_

Model left: \_\_\_\_\_ Serial no.: \_\_\_\_\_

Model right: \_\_\_\_\_ Serial no.: \_\_\_\_\_

## Your individual hearing aid settings

To be filled out by your hearing care professional.

Tinnitus SoundSupport: Limitation on use			
<input type="checkbox"/>	No limitation on use		
	Program	Start-up volume (Tinnitus)	Max volume (Tinnitus)
<input type="checkbox"/>	1	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	2	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	3	Max _____ hours per day	Max _____ hours per day
<input type="checkbox"/>	4	Max _____ hours per day	Max _____ hours per day

## Tinnitus SoundSupport settings overview for your hearing aid

Left			Right	
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Tinnitus SoundSupport	<input type="checkbox"/> Yes	<input type="checkbox"/> No

**A) How to change Tinnitus SoundSupport volume in each ear separately**  
To **increase** or **decrease** the volume (on one hearing aid only), use a short press on the upper or lower part of the push-button repeatedly until you reach your desired level.

**B) How to change Tinnitus SoundSupport volume in both ears simultaneously**

You can use one hearing aid to increase/decrease the sound in both hearing aids. When changing the volume in one hearing aid, the volume on the other hearing aid follows.

To **increase** volume, use a short press on the upper part of the push-button repeatedly.

To **decrease** volume, use a short press on the lower part of the push-button repeatedly.

To be filled out by your hearing care professional.

## General settings overview for your hearing aid





Left			Right	
<input type="checkbox"/> Yes	<input type="checkbox"/> No		<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Short press		Change program	<input type="checkbox"/> Short press	
<input type="checkbox"/> Long press			<input type="checkbox"/> Long press	
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Change volume	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<input type="checkbox"/> Yes	<input type="checkbox"/> No	Mute (miniRITE T only)	<input type="checkbox"/> Yes	<input type="checkbox"/> No
<b>Volume control indicators</b>				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at min/max volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps when changing volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<input type="checkbox"/> On	<input type="checkbox"/> Off	Beeps at start-up volume	<input type="checkbox"/> On	<input type="checkbox"/> Off
<b>Battery indicators</b>				
<input type="checkbox"/> On	<input type="checkbox"/> Off	Low battery warning	<input type="checkbox"/> On	<input type="checkbox"/> Off

To be filled out by your hearing care professional.

## Sound and LED light indicators

Different sounds and light patterns indicate the hearing aid status. The different indicators are listed on the following pages. For light indicators on your charger, see the charger's instructions for use.

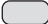
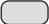
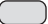

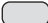




Your hearing care professional can set sound and LED light indicators to match your preferences.

Program	<input type="checkbox"/> Sound	<input type="checkbox"/> LED light*	When to use
1	One tone		
2	Two tones		
3	Three tones		
4	Four tones		

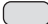

















 Short GREEN blink

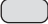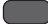\*LED light blinks continuously or is repeated three times with short pauses.

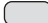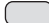Continues on next page

ON/OFF	Sound	LED light	LED light comments
ON	<input type="checkbox"/> Jingle	<input type="checkbox"/>   	Shown once
OFF (miniRITE R only)	<input type="checkbox"/> 4 descending tones	<input type="checkbox"/> 	
Volume	Sound	LED light	
Start-up volume	<input type="checkbox"/> Two beeps	<input type="checkbox"/> 	
Minimum/maximum volume	<input type="checkbox"/> Three beeps	<input type="checkbox"/> 	
Volume up/down	<input type="checkbox"/> One beep	<input type="checkbox"/> 	
Mute		<input type="checkbox"/>  	Continuous or repeated three times

 Long GREEN blink
  Short GREEN blink
  Long ORANGE blink
  Short ORANGE blink

Accessories	<input type="checkbox"/> Sound	<input type="checkbox"/> LED light	LED light comments
TV Adapter 3.0	Two different tones	 	Continuous or repeated three times
EduMic	Two different tones	 	
ConnectClip	Two different tones	 	
Flight mode	<input type="checkbox"/> Sound	<input type="checkbox"/> LED light	
Flight mode activated (miniRITE R only)	Four descending tones + short jingle	   	
Flight mode deactivated (miniRITE R only)	Four descending tones + short jingle	  	Only available when three-time repetition is selected
Flight mode activated (miniRITE T only)	Short jingle	  	Continuous or repeated three times
Flight mode deactivated (miniRITE T only)	Short jingle	 	Only available when three-time repetition is selected



 Long GREEN blink
  Short GREEN blink
  Long ORANGE blink
  Short ORANGE blink

Continues on next page

Warnings	Sound	LED light	LED light comments
Low battery	<input type="checkbox"/> Three alternate tones	<input type="checkbox"/> ● ● ● ● ● ● ● ●	Continuously blinking
	<input type="checkbox"/> Voice Notification		
Battery shut down	Four descending tones		
Microphone service check needed	Eight beeps repeated four times	<input type="checkbox"/> ●●● ●●● ●●● ●●●	Repeated four times
The hearing aid LED light does not turn ON when the hearing aid is placed in the charger <b>(miniRITE R only)</b>		Turned OFF	See the <b>Troubleshooting</b> section

262831AU / 2022.11.12 / v1



Warnings	Sound	LED light	LED light comments
The hearing aid LED light blinks ORANGE when the hearing aid is placed in the charger <b>(miniRITE R only)</b>			Continuously blinking. See the <b>Troubleshooting</b> section.
The hearing aid LED light blinks GREEN when the hearing aid is placed in the charger <b>(miniRITE R only)</b>			Continuously blinking. See the <b>Troubleshooting</b> section.

 Long GREEN blink
  Short GREEN blink
  Long ORANGE blink
  Short ORANGE blink



**oticon**  
life-changing **technology**

262831AU / 2022.11.12 / v1  
Oticon Real Firmware 1.0