



	More 1	More 2	More 3	
Speech Understanding	MoreSound Intelligence™	Level 1	Level 2	Level 3
	- Environment configuration	5 Options	5 Options	3 Options
	- Virtual Outer Ear	3 Configurations	1 Configuration	1 Configuration
	- Spatial Balancer	100%	60%	60%
	- Neural Noise Suppression, Difficult / Easy	10 dB / 4 dB	6 dB / 2 dB	6 dB / 0 dB
	- Sound Enhancer	3 Configurations	2 Configurations	1 Configuration
	MoreSound Amplifier™	•	•	•
	Feedback Prevention	MoreSound Optimizer™ & Feedback shield	MoreSound Optimizer™ & Feedback shield	MoreSound Optimizer™ & Feedback shield
	Spatial Sound™	4 Estimators	2 Estimators	2 Estimators
	Soft Speech Booster	•	•	•
Sound Quality	Frequency lowering	Speech Rescue™	Speech Rescue™	Speech Rescue™
	Clear Dynamics	•	•	-
	Better-Ear Priority	•	•	-
	Fitting Bandwidth*	10 kHz	8 kHz	8 kHz
	Bass Boost (streaming)	•	•	•
	Processing Channels	64	48	48
	Listening Comfort	Transient Noise Management	4 configurations	3 configurations
Wind Noise Management		•	•	•
Personalisation & Optimising Fitting	Fitting Bands	24	20	18
	Multiple Directionality options	•	•	•
	Adaptation Management	•	•	•
	Fitting Formulas	VAC+, NAL-NL1/ NAL-NL2, DSL 5.0	VAC+, NAL-NL1/ NAL-NL2, DSL 5.0	VAC+, NAL-NL1/ NAL-NL2, DSL 5.0
Connecting to the world	Hands-free communication**	•	•	•
	Direct streaming***	•	•	•
	Oticon ON app & Oticon RemoteCare app	•	•	•
	ConnectClip	•	•	•
	EduMic	•	•	•
	Remote Control 3.0	•	•	•
	TV Adapter 3.0	•	•	•
	Phone Adapter 2.0	•	•	•
	Tinnitus SoundSupport™	•	•	•
	CROS/BiCROS support	•	•	•

* Bandwidth accessible for gain adjustments during fitting
 ** Available for Oticon More from FW 1.3 with selected iPhone models
 *** From iPhone, iPad, iPod touch, and selected Android™ devices

Operating and charging conditions
 Temperature: +5°C to +40°C (41°F to 104°F)
 Relative humidity: 5% to 93%, non-condensing
 Atmospheric pressure: 700 hPa to 1060 hPa

Storage and transportation conditions
 Temperature and humidity should not exceed the below limits for extended periods during transportation and storage.

Transport
 Temperature: -20°C to +60°C (-4°F to 140°F)
 Relative humidity: 5% to 93%, non-condensing
 Atmospheric pressure: 700 hPa to 1060 hPa

Storage
 Temperature: -20°C to +30°C (-4°F to 86°F)
 Relative humidity: 5% to 93%, non-condensing
 Atmospheric pressure: 700 hPa to 1060 hPa

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

Oticon More™ miniBTE R is small in size and fits most ears. It is powered by a rechargeable lithium-ion battery. The style features telecoil and a single push-button. It is a Made for iPhone hearing aid and compatible with the new Android protocol for Audio Streaming for Hearing Aids (ASHA) – making it possible to stream directly from iPhone, iPad, iPod touch and selected Android™ devices.

MoreSound Intelligence™ creates a more precise and natural representation of individual sounds with clearer and more distinct contrasts providing access to all relevant sounds.

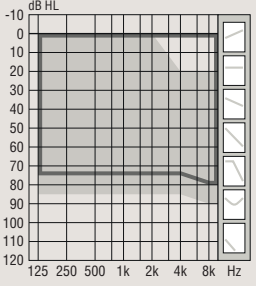

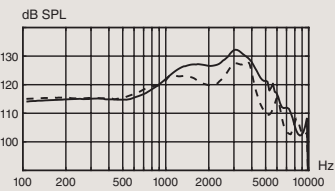
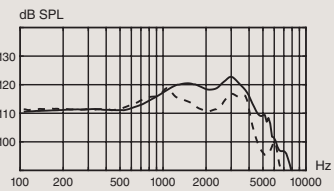
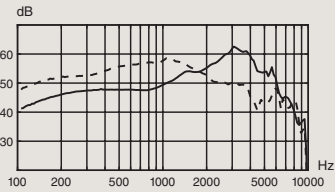
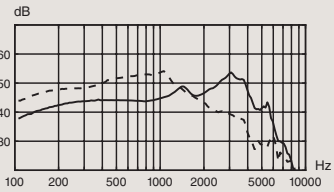
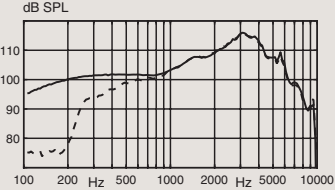
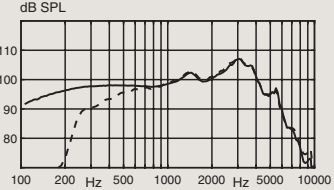
MoreSound Amplifier™ analyses details in sound, and optimally amplifies them for the brain to have access to relevant information.

Oticon More is built on the innovative Polarix™ platform, which uses a Deep Neural Network to rapidly and optimally manage incoming sounds based on individual needs. New features can be added and updates performed wirelessly.



For information on compatibility, please visit www.oticon.global/compatibility



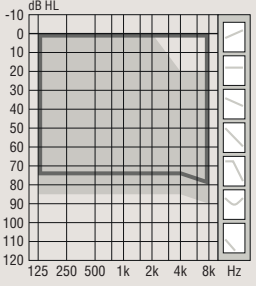

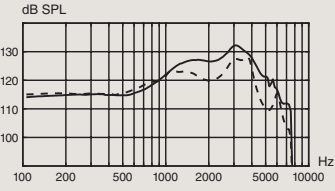
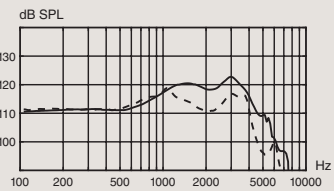
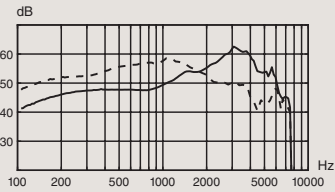
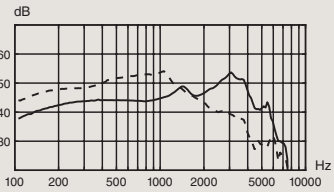
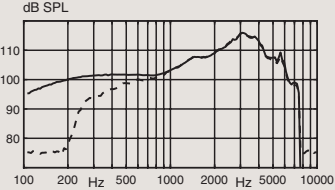
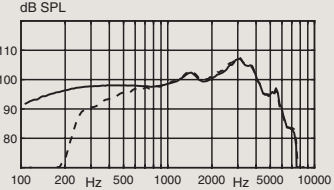
		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 85		OSPL90 	OSPL90 
		Full-on gain 	Full-on gain 
		Frequency response 	Frequency response 
	Peak	132 (128 ¹) dB SPL	123 (119 ¹) dB SPL
	1600 Hz	127 (122 ¹) dB SPL	120 (114 ¹) dB SPL
	HFA-OSPL90	126 (122 ¹) dB SPL	119 (115 ¹) dB SPL
	Peak	63 (59 ¹) dB	54 (54 ¹) dB
	1600 Hz	54 (55 ¹) dB	47 (46 ¹) dB
	HFA-FOG	54 (54 ¹) dB	47 (47 ¹) dB
	Reference test gain	47 dB	41 dB
	Frequency range	100-9500 Hz	100-7300 Hz
	Telecoil output	1 mA/m field (1600 Hz) 10 mA/m field (1600 Hz)	- -
	HFA SPLITS L/R	-	99/99 dB SPL
	Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 4 % < 3 % < 2 %
	Equivalent input noise level	Omni Dir	19 dB SPL 30 dB SPL
	Battery	Lithium-ion	Lithium-ion
	Expected operating time, hours ³	24	

1) For instruments fitted with Corda miniFit Power.

2) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB.

This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.

3) Expected operating time for rechargeable battery depends on use pattern, active feature set, hearing loss, sound environment, battery age and use of wireless accessories.

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 85		OSPL90 	OSPL90 
		Full-on gain  — Standard tube - - - Thin tube (size 1.3)	Full-on gain  — Standard tube - - - Thin tube (size 1.3)
		Frequency response  — Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m	Frequency response  — Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m
OSPL90	Peak	132 (128 ¹) dB SPL	123 (119 ¹) dB SPL
	1600 Hz	127 (122 ¹) dB SPL	120 (114 ¹) dB SPL
	HFA-OSPL90	126 (122 ¹) dB SPL	119 (115 ¹) dB SPL
Full-on gain ²	Peak	63 (59 ¹) dB	54 (54 ¹) dB
	1600 Hz	54 (55 ¹) dB	47 (46 ¹) dB
	HFA-FOG	54 (54 ¹) dB	47 (47 ¹) dB
Reference test gain		47 dB	41 dB
Frequency range		100-7500 Hz	100-7300 Hz
Telecoil output	1 mA/m field (1600 Hz)	85 dB SPL	-
	10 mA/m field (1600 Hz)	105 dB SPL	-
	HFA SPLITS L/R	-	99/99 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz	< 4 %	< 4 %
	800 Hz	< 4 %	< 3 %
	1600 Hz	< 2 %	< 2 %
Equivalent input noise level	Omni	19 dB SPL	17 dB SPL
	Dir	30 dB SPL	32 dB SPL
Battery		Lithium-ion	Lithium-ion
Expected operating time, hours ³		24	

1) For instruments fitted with Corda miniFit Power.

2) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB.

This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.

3) Expected operating time for rechargeable battery depends on use pattern, active feature set, hearing loss, sound environment, battery age and use of wireless accessories.

Headquarters
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark



SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark