

# PC850

Portable PC-based automatic audiometer

The **Amplivox PC850** is a portable automatic audiometer that integrates with a PC to provide a comprehensive and pro-active audiometric facility.

## Comprehensive

Three extensive tests including the Bekesy test and an on-screen audiometric questionnaire facilitate a completely paper-less audiometric function. All tests can be reviewed on-screen and subsequently printed or electronically exported if required.

## Pro-active


The **Amplivox PC850** has a user definable recall facility and extensive data analysis features that enable 'at risk' individuals and user defined groups to be identified and monitored.

## Portability

With a unit weight of just 710g and USB connectivity, the **Amplivox PC850** can be transferred from location to location with ease.

The **Amplivox PC850** is a cost effective solution for the provision, storage and management analysis of audiometric test results.

- Pro-active PC-based automatic audiometer
- Bekesy test function
- Test results analysis
- User definable recall facility
- Identification of 'at risk' personnel
- In-depth management analysis
- Patient communication facility
- USB connectivity

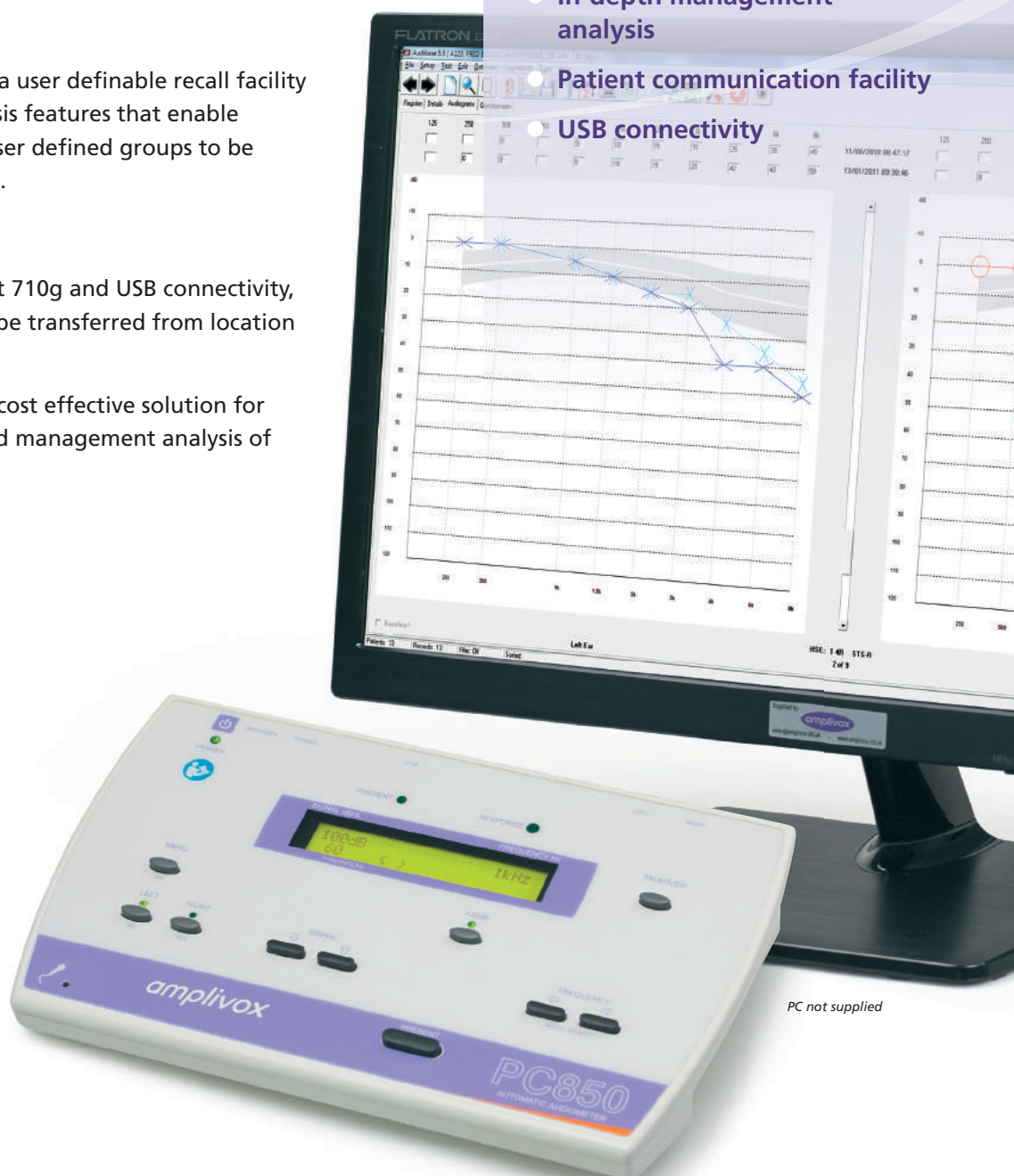

amplivox

specialists in healthcare products

# PC850

Portable PC-based automatic audiometer



## Technical specifications

Output/input data	
Frequency range (Hz):	250, 500, 1K, 1K5, 2K, 3K, 4K, 6K, 8K
Frequency accuracy:	<1%
Distortion:	<2%
Output level range:	-10dBHL to 100dBHL
Output level accuracy:	±3dB
Test modes:	Bekesy Computer (BS6655) Manual
Output level step size:	Bekesy; 2.5dB Computer; 5dB Manual; 5dB
Threshold level:	Bekesy; 1dB Computer; 5dB Manual; 5dB

## System requirements

Operating system:	Stand alone PC with Windows XP or higher
Disc drives:	CD and hard disc drive
Memory:	Internal memory and available disc space as required by the PC operating system
Software:	The programme will be supplied on a CD
Interface:	USB

## Physical data

Display:	2 lines of 24 characters
Power:	100-240V ac 50/60Hz
Dimensions (mm):	270 x 175 x 68 high
Weight:	710 grams

## Standards

The Amplivox PC850 conforms to the relevant clauses of the following standards.

Safety:	IEC 60601-1
EMC:	IEC 60601-1-2
Audiometric:	IEC 60645-1
CE mark:	Complies to EU Medical Device Directive

## Equipment

### Standard equipment

- Audiometric headset
- Patient response switch
- Operation manual
- Mains adaptor
- Quick reference operating guide
- USB cable
- Custom-designed carrying case



### Carrying case

For domiciliary use, the **Amplivox PC850** is supplied with a carrying case, designed to accommodate the audiometer and all accessories.

### Optional equipment

- Audiocups (noise reducing headset enclosures)

**amplivox**

Amplivox Ltd, 6 Oasis Park, Eynsham,  
Oxfordshire OX29 4TP United Kingdom  
Tel: +44 (0) 1865 880846 • Fax: +44 (0) 1865 880426  
E-mail: sales@amplivox.ltd.uk

[www.amplivox.ltd.uk](http://www.amplivox.ltd.uk)