

INTRODUCING
OTICON INO



ENJOY
THE ESSENTIALS



Oticon Ino takes the entry-level hearing solution to the next level. With Ino your clients can experience the benefits of some of the very latest advances in hearing care technology: fast and accurate sound processing, easy connectivity to phones and media devices, and robust feedback suppression - all in a discreet and affordable solution. We have advanced the basic hearing solution, so your clients can interact more naturally, communicate spontaneously and have the confidence to participate actively in all lifestyles, whatever their budget.

ENJOY **THE ESSENTIALS**


ESSENTIAL DETAILS FOR BETTER INTERACTION

Oticon Ino will bring your clients closer to friends and family, and connect them to the world around them. The distinct sound quality provided by our latest RISE 2 platform lets your clients enjoy the details and richness of the sounds around them. Access to

details is essential to communication and relaxed interaction - the listening experience becomes less strenuous when the nuance of voices come through more clearly.

Oticon Ino is built on Oticon's most powerful processing platform to date: **RISE 2**, enabling:

- Double processing speed compared to RISE
- Great sound quality
- Detailed sound picture
- 8 kHz bandwidth
- Dynamic Feedback Cancellation

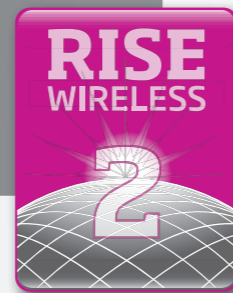



ESSENTIAL CONFIDENCE NOW AND FOR THE FUTURE

Oticon Ino is a sound investment for now and for the future. It is available in a range of small and lightweight designs so your clients can choose precisely the shape and style that suits their needs - and feel comfortable and confident whenever they are using it.

Ideal for many of your clients, Oticon Ino is exceptionally easy to fit, and also features in-situ audiometry which lets you use the clients' own moulds and hearing instruments when measuring the hearing threshold. Should your client's hearing needs change, the speakers can be replaced to obtain higher outputs to cater for the new requirements.



Oticon Ino is also available in one of our most discreet styles to date, **miniRITE**


- Tiny and lightweight
- Push button for easy control
- Fits comfortably and discreetly behind the ear
- Advanced wireless capabilities
- Fits up to 110dB HL
- Nano-coated for endurance

Oticon Ino is remarkably easy to fit and to adapt:

- Broad range of styles, including the ultra-discreet miniRITE
- Can be fitted to a wide range of users and hearing losses
- Flexible in-situ audiometry
- Wireless fitting option
- Flexible RITE speaker system offering 3 power levels

**ESSENTIAL ACCESS
TO COMMUNICATION AND
ENTERTAINMENT**

Oticon ConnectLine is the complete connectivity system, which empowers people in more situations than ever before. Barriers to natural communication disappear, as the user intuitively participates, just by touching a button on the ConnectLine Streamer.

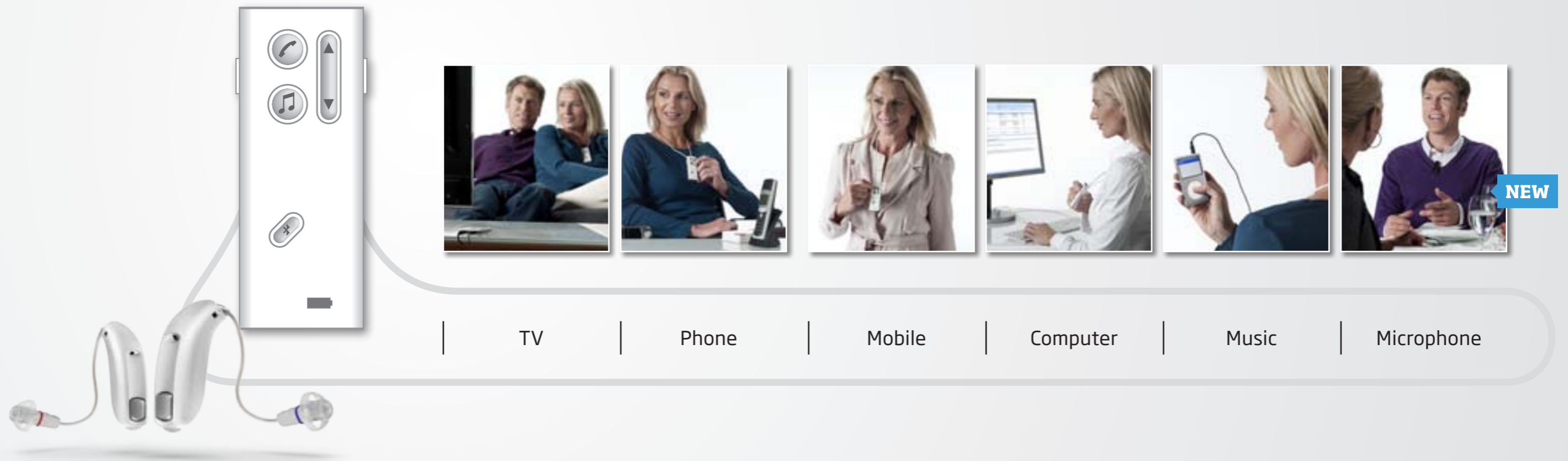

ConnectLine TV

With ConnectLine TV, the user is provided with a TV sound solution that is better than with just hearing instruments alone. For the user, this means high sound quality and easy, personalized volume control. For the rest of the family, it means watching TV at their choice of volume level.

ConnectLine Phone

Stress-free conversations are now possible thanks to ConnectLine Phone. The Phone adapter wirelessly connects to the Streamer, so that users can take calls remotely with the mere touch of a button. The Streamer has a built-in microphone - in effect turning a pair of hearing instruments into a wireless handsfree solution for the home phone.

Mobile

The Oticon Streamer can connect directly to any Bluetooth™ mobile phone. The built-in microphone on the Streamer allows users to answer a call handsfree - without having to touch, or even find, the phone.

Computer

More and more people are using PCs for internet telephony (Skype, FaceTime, Messenger, etc.), video conferencing, watching films and much more. The Streamer can connect to the PC in two ways: via Bluetooth or using the jack cable supplied.

Music

Connect Streamer to radios, PCs, MP3 players and other audio devices either via Bluetooth™ or with the cable supplied - and change the hearing instruments into headphones.

ConnectLine Microphone

The ConnectLine microphone is a discreet microphone that connects wirelessly to the Streamer and lets users enjoy clear and private one-to-one communication in challenging listening environments.

ESSENTIAL OVERVIEW

miniRITE/RITE



miniRITE

RITE

STANDARD / MEDIUM



Open Dome



Plus Dome



Micro Mould

POWER




Power Dome



Power Micro Mould



Power Mould

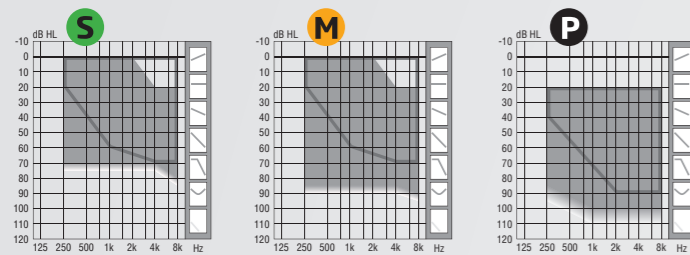

STANDARD

MEDIUM

POWER

■ Custom Mould - Power Dome

□ Open Dome - Plus Dome

CUSTOM



ITE*



ITC Power (Omni or Dir)**



ITC*




CIC/MIC Power**



CIC/MIC**

* Also available as a non-wireless style
** Not wireless

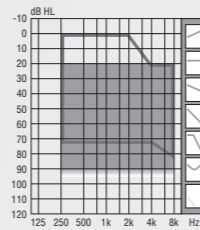

■ ITC Power

□ CIC/MIC


BTE 13



Corda² thin tube



Hook

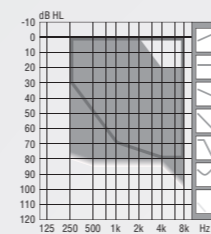


■ Earmould

□ Corda2

BTE POWER



Hook

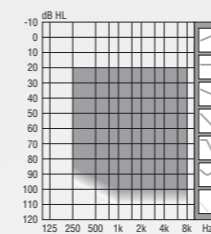

Main Features

- Binaural Coordination
- Dynamic Feedback Cancellation (DFC2)
- AI Essential
- Adaptive Directionality
- Automatic Directionality (Surround and Split)
- Noise Management

Ino Pro

-
-
-
-
-
-

Ino

-
-
-
-
-
-



Fitting options (Genie 2011.1)


nEARcom enabled	•	•
In-situ audiometry	•	•
NAL-NL1, NAL-NL2, DSL v5.0a m[i/o]	•	•
Identities	1	1
Fitting Bands	6	4
Adaptation Manager	Auto	Manual

More

Bandwidth	8 khz	8 khz
3 RITE receivers (S, M and P)	•	•
RISE 2 platform	•	•

Connectivity

ConnectLine enabled	•	
DAI and FM	•	•

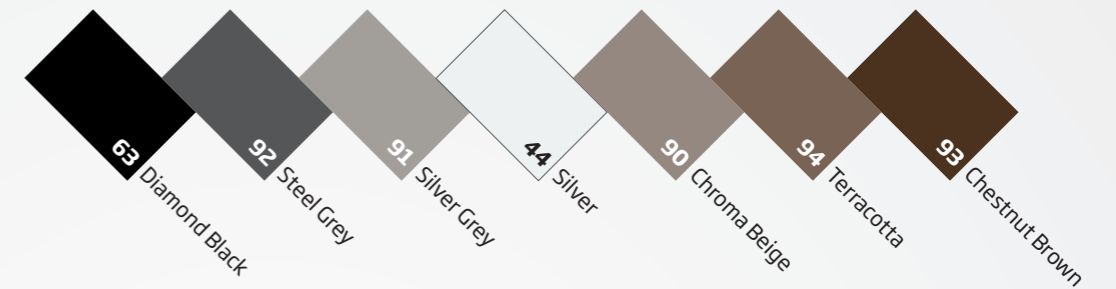

People First

People First is our promise
to empower people
to communicate freely,
interact naturally and
participate actively