

Otoport

TE SCREENER

Ideal For Universal Newborn Hearing
Screening Programmes



Australian Distributor



SONIC
EQUIPMENT
One stop shop for all your audiological needs

1800 639 263
info@soniceq.com
soniceq.com





Otoport

TE SCREENER

Ideal For Universal Newborn Hearing Screening Programmes

Otoport Screener is a fast TE Newborn and Infant Screener, which integrates with Universal Newborn Hearing Screening Programmes.



FULLY CUSTOMISABLE
3 Protocols, 16 Parameters



QUALITY CONTROL
Dedicated QA Area



FREQUENCY RANGE
1 to 6 kHz bands, $\frac{1}{2}$ oct





OVERVIEW


- Light-weight and portable with easy-to-use keypad
- Automated self-test system when powering on Otoport
- Intuitive test progress graphic displaying signal to noise progress
- Long battery life - up to 250 test cycles or one full working week
- On-screen notification of poor probe fit, high noise level and ringing stimulus
- Each test is saved with the actual physiological response
- Comprehensive probe fit and noise level indicators
- Create multiple User logins with password protection with Admin control if required
- Fast, noise-resistant mode
- On-board probe test facility and Probe ID tracking
- Transfer data wirelessly to the Otolink suite
- Print wirelessly to the Otodynamics optional mini printer

TE OAE	
Test Configurability	3 Protocols, 16 Parameters
Frequency Range	1, 1.5, 3, 4, 5, 6kHz
Freq. Resolution	1/2 octave bandwidth
Stimulus Levels	60 to 90 ($\Delta 1\text{db}$) (dB SPL _{pe})
Pass Criteria	7 Parameters
Pass SNR	1 to 14dB (Fixed)




FAST TEOAE SCREENING


For ultimate test documentation, the Otoport Screener is the ideal for Universal Hearing Screening Programmes providing TEOAE Quickscreen with customisable protocols, configurable individual pass criteria, dedicated programme QA area and QA history storage. Plus, searchable 5000 test and patient database with extensive biographic data fields.

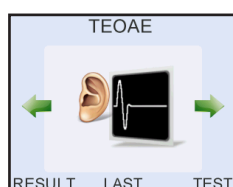


Easy 'press-and-go' operation with choice of 2 key press test start.

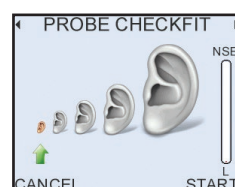


In-ear calibration ensuring accurate results with customisable stimulus levels 60-90dB 1dB steps.

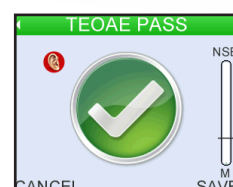


Fast and clear Pass/Refer result. Choice of overall or frequency band displays.

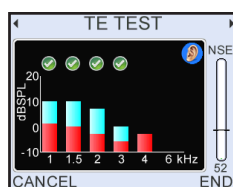

High resolution spectrum display up to 6kHz.

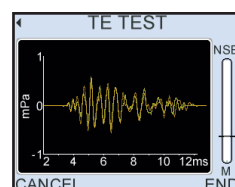

Waveform display, dB SPL, SNR and Repro.

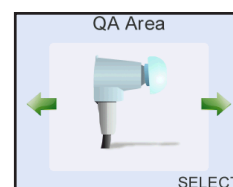

Dedicated QA test area On board probe, cavity, occlusion and real ear test QA facility with Probe ID tracking & history view.



EASY DATA ENTRY & RECALL

WORK LIST 3/3	
NAME	ID/NOTES
Roberts	OGPM1V02
Jones	OGPM1V03
Thomson	OGPM1V04
DELETE	DETAILS TEST

Work List and downloadable secure database. Login access if required.

Grey 1/1	
	TEOAE
	TEOAE PASS
	Any2Bands, 84dB
OAE	7.1 dB
Noise	1.9 dB
24 Aug 17	18:55 ADN
BACK	PCLOAD SELECT

Extensive 5000 test 1000 name record database, searchable by name, ID and date, tester, DoB.

ENTER DETAILS	
ID/Notes	OGPM1V01
Name	30 Dec 2011
D.O.B	31 Jan 2012
Gender	01 Feb 2013
CANCEL	INSERT

Extensive patient database fields Enter Full Patient demographic data as needed.

DATA SECURITY

USERS 1/1	
USER	PASSWORD
Admin	A
BACK	EDIT

Create multiple User logins with password protection with Admin control if required.

USERS	
1 Add New User	
2 View Users	
3 Settings	
BACK	SELECT

Extensive protocol control with login and Administrator functions.






OTOLINK

PC Data Management Suite

Manage your data easily and efficiently with the Otolink PC patient Data Management Suite. Download test results wirelessly or via USB from the Otoport or when docked in the docking station*. Test details can be quickly reviewed, printed and archived. Printing reports directly to PDF file is available. Data can be exported, for example, to OZ, HiTrack, GDT, Noah** or as CSV file. Patient Work Lists can be uploaded to the Otoport, assisting with managing workloads in busy environments. Otolink also provides access to firmware and application updates from Otodynamics via the internet. Compatible with Windows 10, 8 and 7.

* With docking station option ** With Noah option

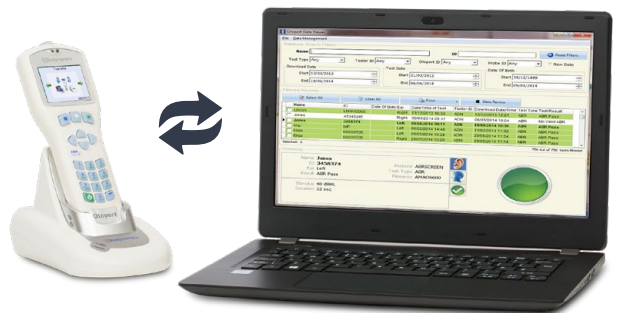

Ez-SCREEN

PC Software

USB transfer to PC for Ez-Screen software analysis and tracked integration into Ez-Screen database, exporting to OZ, HiTrack, ASCII, XML, archive. Compatible with Windows 10, 8 and 7.



WHAT'S INCLUDED

with your Otoport

- Charging cradle
- OAE probe with cable clip
- Sample ear tip pack
- Probe test cavity
- Sample infection control sleeve
- PC Otolink & Ez-Screen software
- USB download cable
- User manuals
- Carry case



2 YEARS
INSTRUMENT
WARRANTY

For further information, please contact us:

www.otodynamics.com
info@otodynamics.com





ADDITIONAL options for your Otoport

OAE + **ABR**

Add AABR with the easy-to-fit sleeve

At any time add ABR screening to your Otoport with the easy-to-fit ABR sleeve. Ideal for Newborn Screening. Configurable dBnHL stimulus level. Choice of pass criteria. Automated start for speed. Clear pass/refer indication. Real time ABR waveform view. Work List and User login option for large screening programmes.

(See Otoport ABR brochure.)



ADD MINI PRINTER

Compact, fast print time, numerical and band pass data. Battery powered. Choice of wired*/ wireless printing via the Otoport.

*Wired printing not compatible
with ABR sleeve



ADD DOCKING STATION

The docking station replaces the Otoport charging cradle with full connectivity to the OAE Otoport, PC, charger and printer with automatic print-on-connect mode.

(Not compatible with ABR sleeve.)



Noah 4

ADD NOAH

Our Noah option makes Otolink PC Data Management Suite fully compatible with Noah 4 based data systems.



Australian Distributor

SONIC
EQUIPMENT
One stop shop for all your audiological needs

1800 639 263
info@soniceq.com
soniceq.com

Otodynamics Ltd.

30-38 Beaconsfield Road Hatfield Herts AL10 8BB UK
31-+44 (0)1707 267540 E. info@otodynamics.com

www.otodynamics.com

Ref: F-008-11ROW iss002